ENVIRON[®]

Focus Care Skin Tech+

Electro-Sonic

DF Mobile Skincare Technology

INSTRUCTION MANUAL

Please read and retain this booklet for future reference and warranty registration.

ENGLISH	3
FRANÇAIS	12
DEUTSCH	21
NEDERLANDS	30
ESPAÑOL	39
PORTUGUÊS	48
ITALIANO	57
POLSKI	66
ΕΛΛΗΝΙΚΑ	75
NORSK	84
SUOMI	93
SVENSKA	102
EESTI KEEL	111
DANSK	120

Instruction Manual

Creating a Beautiful Future with the DF Mobile

THE ULTIMATE IN SKINCARE TECHNOLOGY

Congratulations on purchasing the revolutionary Electro-Sonic DF Mobile and taking your first steps in creating your future beautiful skin.

CONTENTS

DF MOBILE DESCRIPTION	4
CONTRAINDICATIONS	4
POTENTIAL SIDE EFFECTS	5
WARNINGS	5
PRECAUTIONS	5
BEFORE THE FIRST USE	6
DIRECTIONS FOR USE	7
DF MOBILE CARE	9
TECHNICAL DETAILS	11

The Ionzyme® DF Mobile Machine is not intended for use in the diagnosis of disease or other conditions, or in the cure, mitigation, treatment or prevention of disease, and is not intended to affect the form, structure or any function of the living tissue and therefore is not classified as a device under the FD&C Act. 179.


Focus Care Skin Tech+® Electro-Sonic DF Mobile Skincare Technology
Manufactured for: Environ Skin Care (Pty) Ltd
14 Jan Smuts Road, Beaconvale, Parow, Cape Town, 7500, South Africa
Made in South Africa

www.environskincare.com

DF MOBILE DESCRIPTION

The DF Mobile uses low frequency Sonophoresis and Iontophoresis and is intended to temporarily disrupt the Stratum Corneum to improve the penetration of skincare for cosmetic applications. It is intended to be used by a layperson in the home environment. Do not use the DF Mobile for any other intended purpose.

By temporarily disrupting the Stratum Disjunctum (looser and superficial layer of the Stratum Corneum) the barrier function created by the Stratum Corneum is reduced (not removed) therefore helping cosmetic skincare products to be absorbed into the Stratum Compactum (lower layers of the Stratum Corneum).

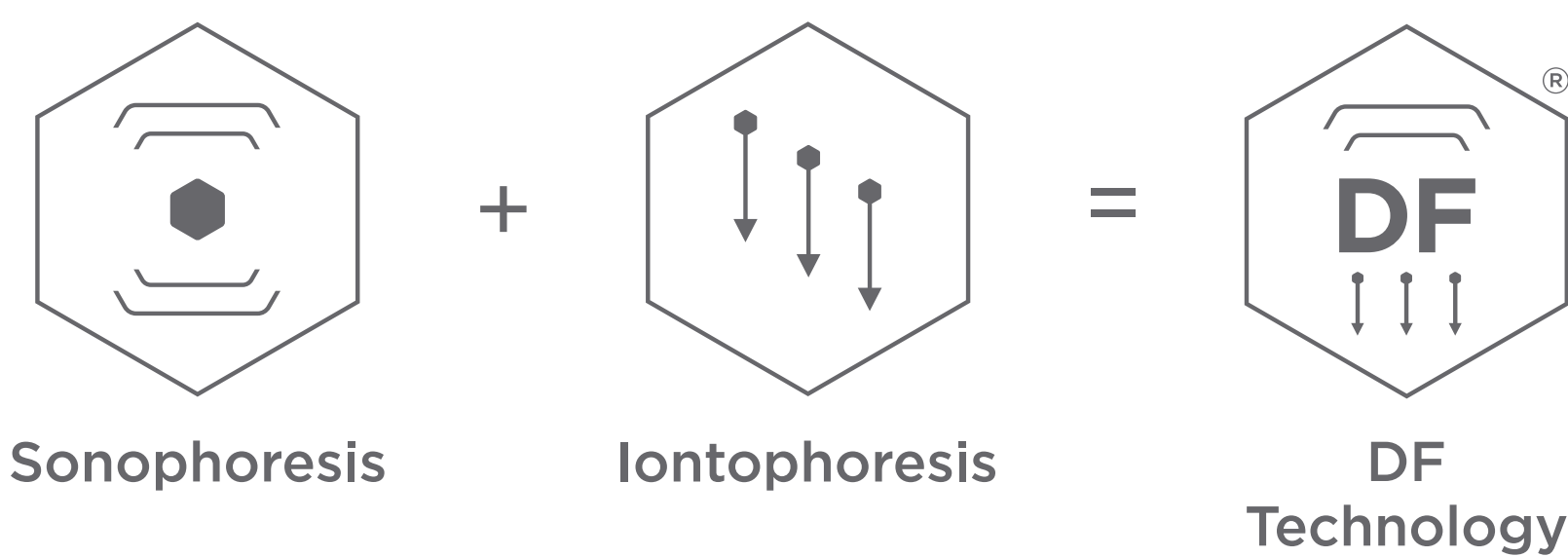Low Frequency Sonophoresis uses specific sound waves while Pulsed Iontophoresis uses a low level electrical current. This dual electro-sonic technology assists in delivering essential nutrients to where skin needs them most.

Combined with Environ's cosmetic products, the appearance of specific skin concerns can be targeted using the DF Mobile, resulting in the look of a noticeably improved complexion, texture and radiance.

To keep the skin looking healthy and more beautiful, select the appropriate programme before application:

- Programme 1** – for enhanced effectiveness and delivery of skincare
- Programme 2** – for enhanced cleansing
- Programme 3** – for normal effectiveness and delivery of skincare with gentle sonic massage
- Programme 4** – for minimum enhanced effectiveness and delivery
- Programme 5** – for mild cleansing

CONTRAINDICATIONS

This technology is contraindicated for the following conditions:

1. Pacemakers and electronic implants of all kinds
2. Pregnant or lactating women
3. Epilepsy
4. Open cuts & abrasions in the treatment area
5. Cold sores (Herpes simplex type I) and ulcers
6. Metallic implants of the face
7. Botulinum Toxin treatment and Dermal Fillers: A two week waiting period is recommended to avoid affecting their placement

8. Sunburnt skin
9. Sufferers of angina/arrythmia
10. Haemophilia
11. Patients who have a history of stroke
12. Cancer patients or those undergoing radiotherapy

POTENTIAL SIDE EFFECTS

1. Erythema – superficial reddening of the skin may be observed which subsides soon after the treatment has been completed.
2. No other side effects have been reported.

WARNINGS



1. Always unplug the charger cable before use. The DF Mobile will not work while connected to the power source.
2. Do not immerse the DF Mobile in water, other liquids or expose to moisture.
3. Always unplug the DF Mobile if there has been any danger of water or other liquid spillage on to it.
4. Ensure that the conduction cloth is only slightly damp to avoid damaging the DF Mobile and to ensure effective and safe conduction of the electrical current.
5. Only use the supplied charger to charge the DF Mobile.
6. Use of the DF Mobile is only intended for use as described in this instructional manual.
7. Never charge the DF Mobile if it has a damaged cable or plug, if it is not working properly, or if it has been dropped or damaged.
8. Keep the charger cable away from heated surfaces.
9. Unplug the DF Mobile from the wall outlet if you do not use it for a long time. Disconnect at plug point. Never pull on the cable itself.
10. In lightening-prone areas, unplug the DF Mobile from the wall outlet during an electric storm.
11. Stop the treatment immediately should irregularity or abnormalities of the DF Mobile or the skin occur.
12. Do not tamper with the DF Mobile or attempt to modify it in anyway.

PRECAUTIONS


Remove all jewellery before using the DF Mobile.

BEFORE THE FIRST USE


Unpack the contents of the carton and familiarise yourself with the DF Mobile.

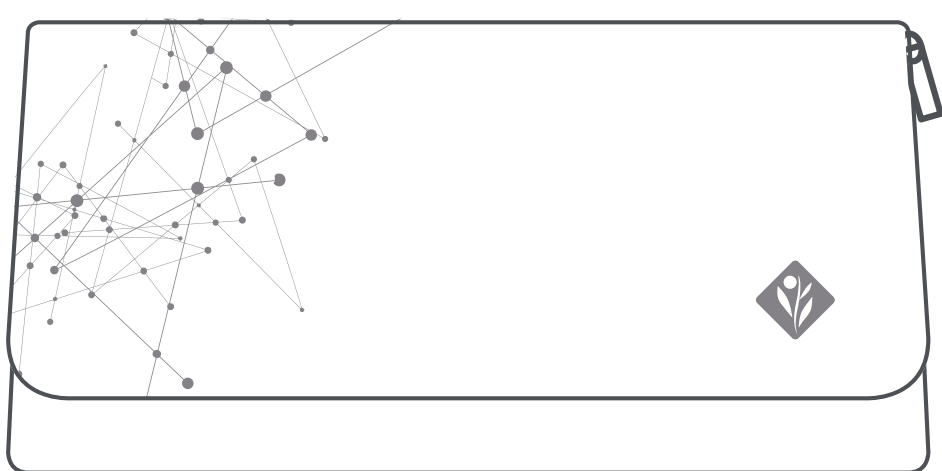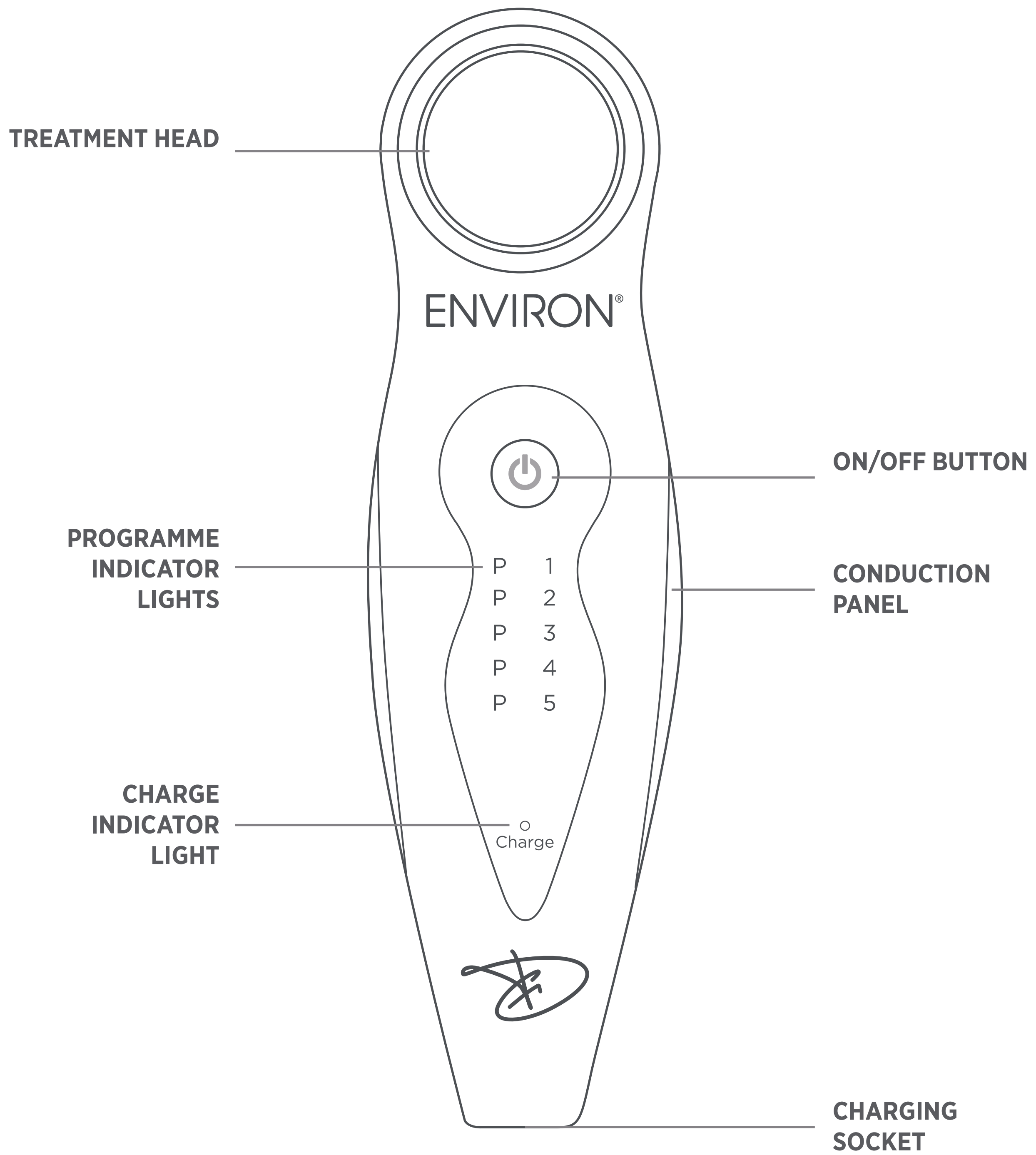
TRAVEL BAG

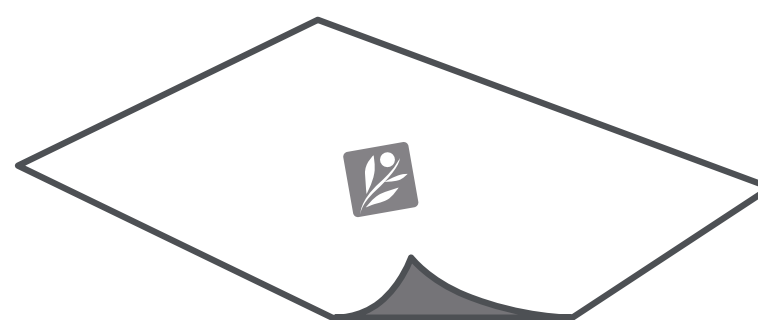
CONDUCTION CLOTHS

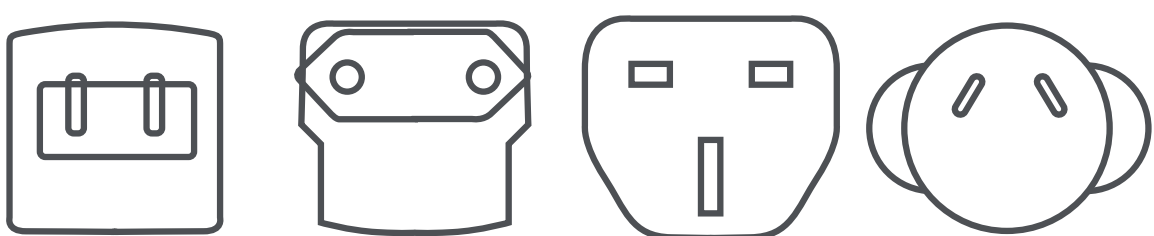INTERNATIONAL PLUGS

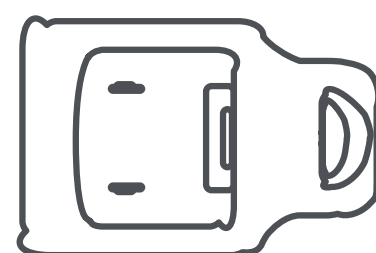CHARGER

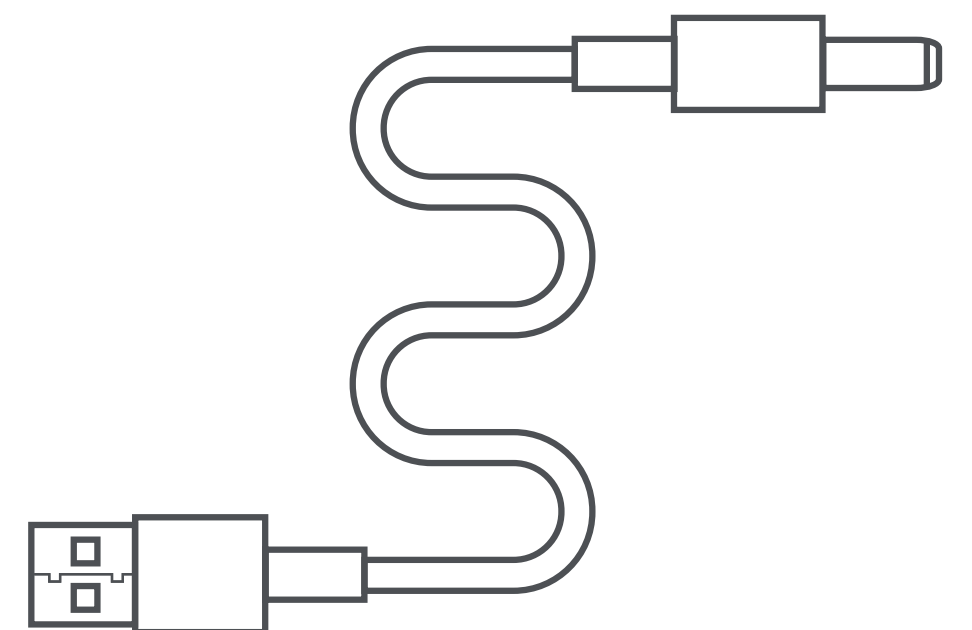USB CHARGING CABLE

HOW TO CHARGE/RECHARGE YOUR DF MOBILE

1. Connect the **CHARGER** into the required **INTERNATIONAL PLUG**.
2. Connect the **USB CHARGING CABLE** to the **CHARGER**.
3. Plug the **CHARGING CABLE** into the **CHARGING SOCKET**. Plug the **CHARGER** into wall socket and switch on at the wall socket.
4. The **CHARGE INDICATOR LIGHT** will illuminate orange and remains orange while being charged. Maximum charge time will be 45 – 60 minutes. When the **CHARGE INDICATOR LIGHT** changes to green, it indicates that the battery is fully charged.

Note: The DF Mobile cannot be used while being charged and the orange light is on – this is a safety feature. Charging the DF Mobile when it has not been fully discharged is safe.

IMPORTANT

1. When the battery is depleted the DF Mobile will not work and no lights will illuminate when pressing the ON/OFF button.
2. Recharge your DF Mobile using your charger as described above.
3. There are no over-charging restrictions. You can leave the charger connected to the DF Mobile when not in use.

DIRECTIONS FOR USE

Once your DF Mobile is fully charged you are ready to begin a treatment.

1. Remove all jewellery before use.
2. Clean the skin thoroughly.
3. Apply a water-based gel to the skin which will ensure good conduction of both electricity as well as sound to the area to be treated. *Selecting the appropriate gel: Ask your Environ trained Skincare Professional to recommend Environ Skin Care product(s) best suited to your specific skin concern.*
4. Dampen the **CONDUCTION CLOTH** and wrap it around the back of the DF Mobile, ensuring that the cloth covers both metal panels. The palm of the hand should make good contact with both panels via the cloth. The cloth must be slightly damp, not wet, to prevent water-damage to it.
5. Apply the same gel to the **TREATMENT HEAD** of the DF Mobile.
6. Press the **ON/OFF BUTTON** for one second to switch the DF Mobile ON.

- Place the flat part of the **TREATMENT HEAD** flush onto the selected area of the skin and wait for a “beep” sound to indicate the start of the treatment.

Hold it in a static position until the next beep sound is heard (15 seconds), before moving it slightly, always slightly overlapping the previous position.


Repeat this until the defined area has been covered. (See diagrams and list below for areas that may be treated) One may treat either small areas or even the whole face. It is unlikely that any individual will use the DF Mobile for longer than 15 to 20 minutes at most. Typically, most people report that they use it for about 5 to 10 minutes in total.

- Repeat the treatments twice a week in the beginning and slowly increase frequency until the product can be used once daily.


WEEK 1	Apply every third day
WEEK 2	Apply every third day
WEEK 3	Apply every second day
WEEK 4	Apply every second day
WEEK 5	Apply daily

TREATMENT AREAS


- Frown lines on forehead and/or between eyebrows
- Smile lines and fine wrinkles around the eyes
- Upper lip lines
- Hyper-pigmented marks (sun damage)
- Minor scars (not indicated for severe scarring/keloid scarring)
- Double chin (submental area)

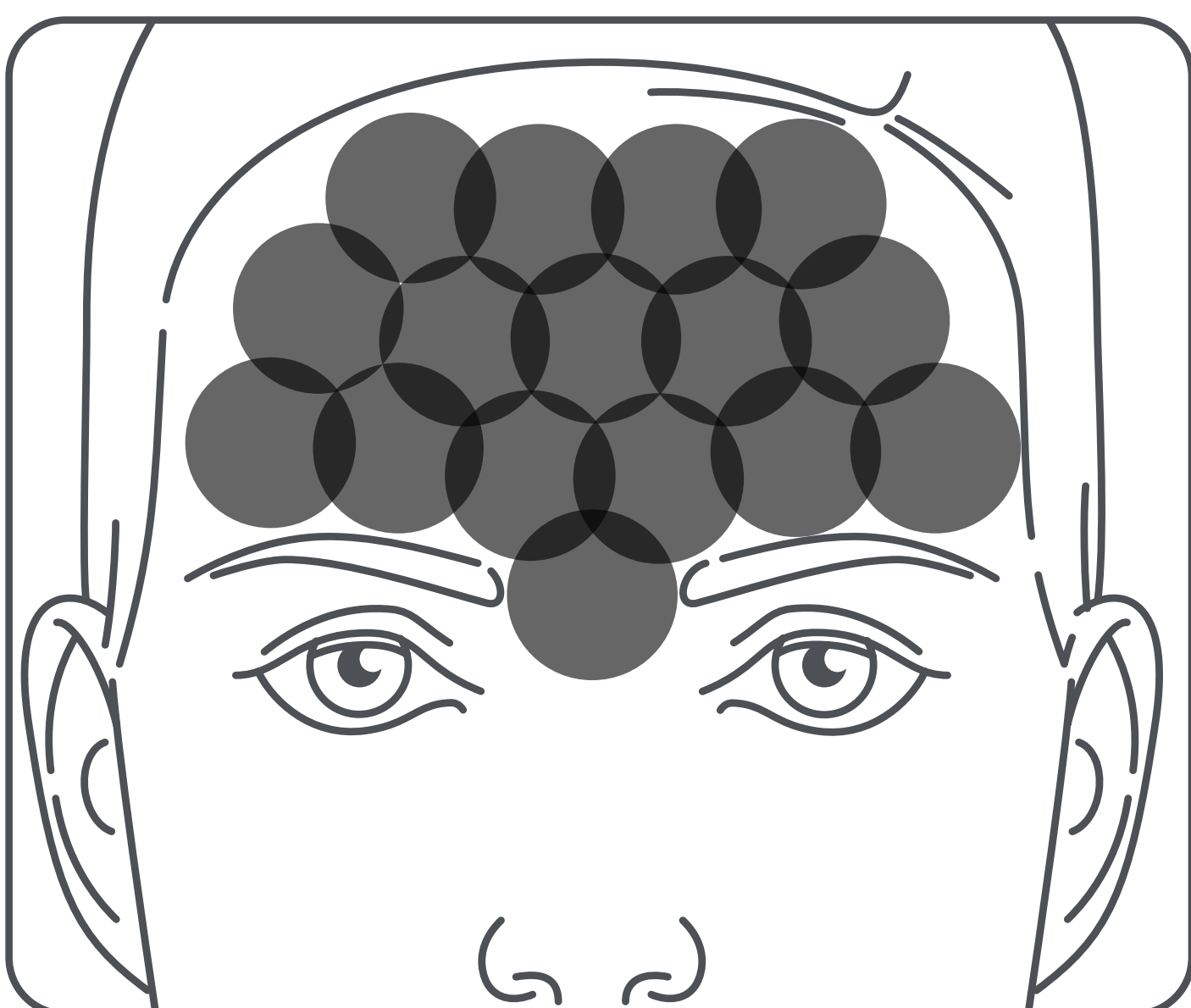
Frown lines on forehead and/or between eyebrows

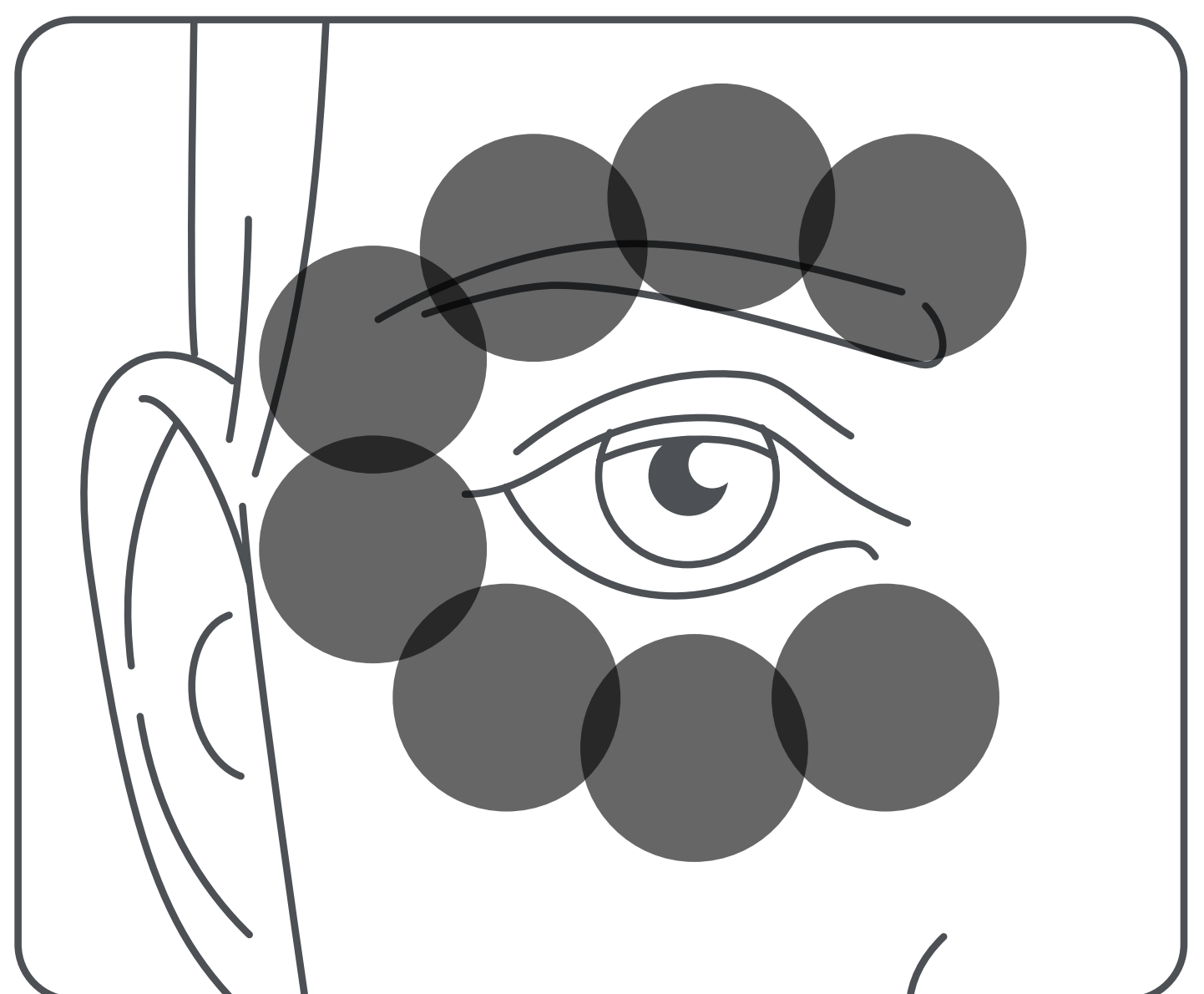Smile lines and fine wrinkles around the eyes

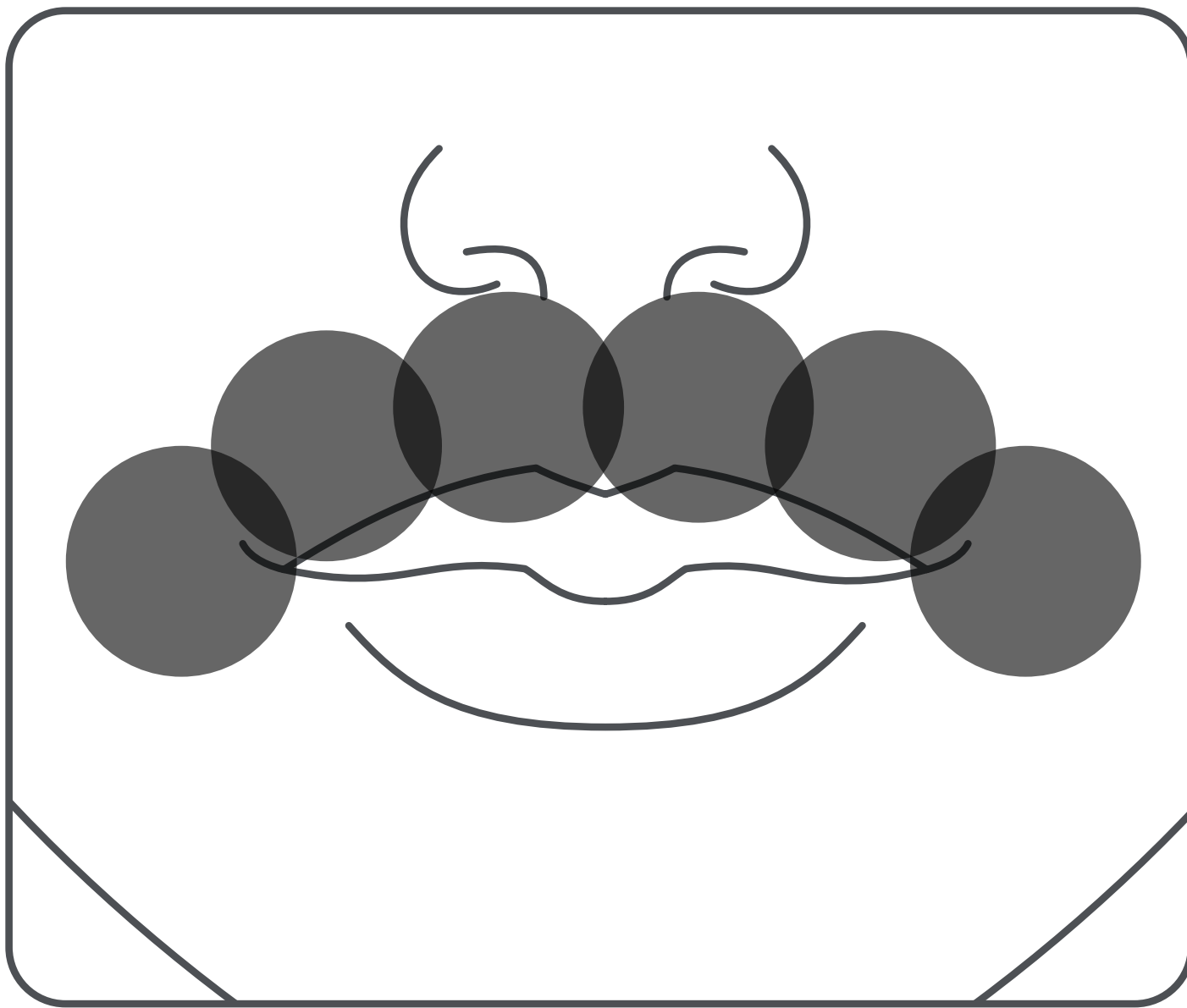Upper lip lines

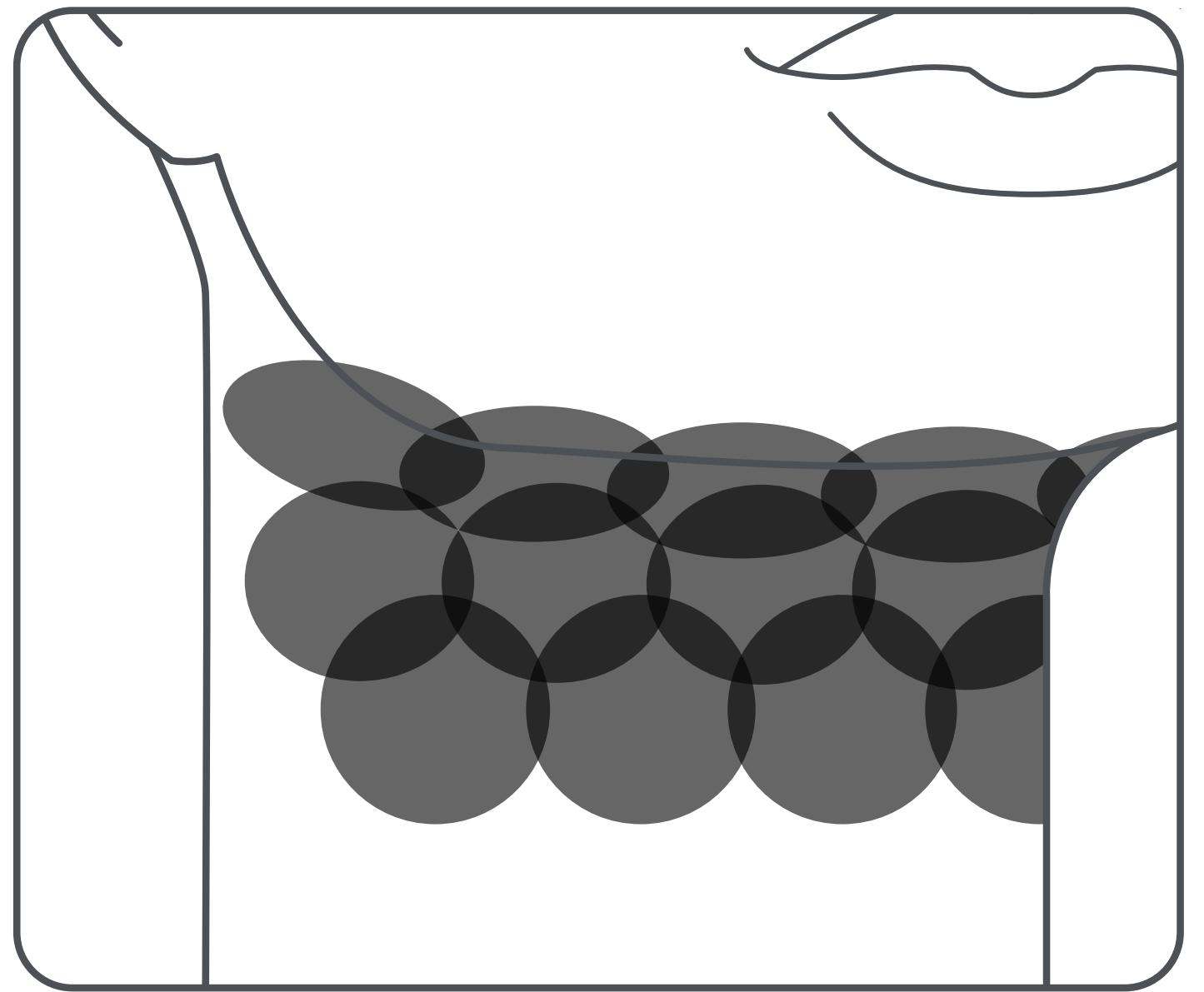Double chin (submental area)

Note: Be sure to always use a slightly damp conductive cloth with your DF Mobile. It is normal to feel a slight tingling sensation in your hand while doing the treatment. This is dependent on the sensitivity and thickness of the skin, rendering greater or less resistance to the electrical part of the treatment. Also, a high pitch oscillating sound may be heard, which is normal.

DF MOBILE CARE

LOOKING AFTER YOUR DF MOBILE

1. Before and after use:
Wipe the DF Mobile with a clean, soft, and slightly damp cloth.
2. Do not immerse the DF Mobile in water.
3. Periodically clean the **TREATMENT HEAD**, by wiping with a soft cloth using rubbing alcohol or 75% alcohol. Do not remove the **TREATMENT HEAD**.
4. Store the DF Mobile safely in the **TRAVEL BAG** when not in use.

BATTERY CARE

1. The DF Mobile automatically switches off after 15 minutes. Each time you start a new session it will switch off automatically after 15 minutes.
2. The DF Mobile should give you at least eight to ten, 15-minute sessions before it requires recharging.
3. If you do not intend on using the DF Mobile for an extended period of time, you will need to do a full recharge of the DF Mobile until the green **CHARGER INDICATOR LIGHT** shines, before storing it.
4. If the DF Mobile will not work and the battery does not charge, check that all the plugs are securely plugged in. The DF Mobile cannot be operated while being charged.
5. To replace the battery, return the machine to any Authorised Environ stockist. They will return it to their Authorised Environ Distributor for replacement.

WARRANTY

Register your DF Mobile to activate the warranty on www.environskincare.com using the unique serial number found on the packaging. Please retain your proof of purchase.

Do not open or attempt to make any modifications to the DF Mobile - any interference will automatically negate the warranty.

SPECIAL PRECAUTIONS

1. Disposal of the DF Mobile - Please return it to a listed Authorised Environ Stockist / Distribution Centre for responsible disposal.
2. If the DF Mobile does not function according to the instructions, please discontinue use immediately and contact a listed Authorised Environ Stockist for assistance. Do not attempt to repair it yourself.

To return the DF Mobile, please follow these instructions:

- Carefully pack it into its original carton, ensuring it is clean.
- Enclose the proof of purchase and your proof of warranty.
(*2-year warranty on the DF Mobile*)
- Provide a reason for it being returned.

TECHNICAL DETAILS

BATTERY:	3,7 V 500 mA
POWER:	0,15 Watt Max
DIMENSIONS:	175mm x 50mm x 48mm
MASS:	110 g
ELECTRICAL CLASSIFICATION:	CLASS 1
MAXIMUM ON TIME:	15 Minutes

CHARGER UNIT

ELECTRICITY:	100-240V AC, 50-60 Hz
POWER:	10 Watt Max
DIMENSIONS:	46 mm x 55 mm x 36 mm w/o adaptor
MASS:	40 g w/o adaptor
ELECTRICAL CLASSIFICATION:	CLASS 1
MAXIMUM ON TIME:	Continuous

AMBIENT TEMP:	10 °C to 40 °C
RELATIVE HUMIDITY RANGE:	30% to 75%, non-condensing
ATMOSPHERIC PRESSURE:	700 hPa to 1060 hPa

STORAGE AND TRANSPORT

AMBIENT TEMPERATURE RANGE:	- 40 °C to + 70 °C
RELATIVE HUMIDITY RANGE:	10 % to 100 %, including condensation
ATMOSPHERIC PRESSURE RANGE:	50 kPa to 106 kPa

SONOPHORESIS

ACOUSTIC WORKING FREQUENCY:	22 kHz $\pm 10\%$
EFFECTIVE RADIATING AREA:	7,07 cm ²

IONTOPHORESIS

OUTPUT DC	pulsed
OUTPUT WAVE FORM	pulsed DC Square Wave with equal ON/ OFF periods of 15 mS, 33Hz $\pm 10\%$
OUTPUT POLARITY	either Positive or Negative (switched)
VOLTAGE 24 V	pulsed $\pm 10\%$ (Maximum 26.4 V)
CURRENT	0.2 mA $\pm 15\%$ (Maximum 0.23 mA)

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