

# HOW TO USE THE DF MOBILE

## *Preparation steps*

- ◇ **STEP 1:** PREPARE THE SKIN
- ◇ **STEP 2:** SELECT THE APPROPRIATE SERUM FOR YOUR SKIN CONCERN

## *Treatment steps*

- ◇ **STEP 3:** SWITCH ON AND TREAT
  - Make sure the DF Mobile is on setting P1.
  - Dampen the conduction cloth and wrap it around the DF Mobile.
  - Apply the selected Environ product to area of concern.
  - Place the flat part of the Treatment Head onto the selected area and wait for a “beep” sound to indicate the start of treatment.
  - Hold the DF Mobile in a static position until the next beep sound is heard (15 seconds) before moving it slightly.
  - For a more intensive session, repeat the process more than once.
  - On completion, massage any excess serum into the skin and follow with a recommended Environ moisturizer and sunscreen.


## *The benefits of using the DF Mobile*

The DF Mobile helps to allow for enhanced effectiveness and delivery of ingredients to where the skin needs them most. The result is an **improvement in the appearance of skin health and wellness** which increases the more the DF Mobile is used.

The DF Mobile offers the advanced skincare user a fantastic tool that is safe and gentle enough to be used every day as part of a comprehensive skincare regime.



# ELECTRO-SONIC DF MOBILE *Skincare Technology*

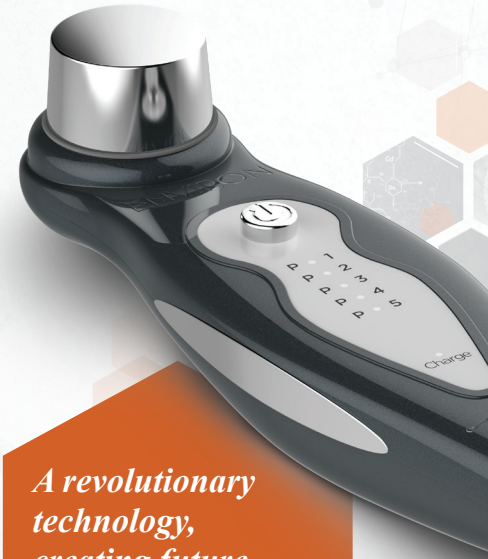


*A revolutionary  
technology,  
creating future  
beautiful-looking  
skin.*

# WHO SHOULD USE THE DF MOBILE

The DF Mobile should be used by anyone who wants to maximize the benefits of professional treatments through the combination of skincare products and a home-use appliance to target the appearance of the signs of the following areas of concerns:

- Frown lines on forehead and/or between eyebrows
- Smile lines and fine wrinkles around eyes
- Upper lip lines
- Hyper-pigmented marks from sun damage
- Minor scars (however, it is not indicated for severe scarring/keloid raised scarring)
- Double chin (submental area)

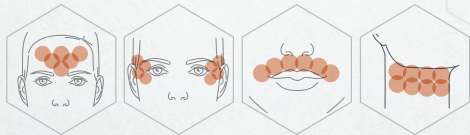

# WHY USE THE DF MOBILE


The DF Mobile delivers the nutrients the skin needs – where skin needs them – making a lasting difference, and giving you the confidence to live life beautifully.

- *Providing nutrients to where the skin needs them most*
- *Delivering ingredients to the skin between professional treatments*
- *Innovative skincare technology in your hands*



# Electro-Sonic DF Mobile

This breakthrough technology makes a lasting difference to the appearance of the skin and allows a user to get better results faster.

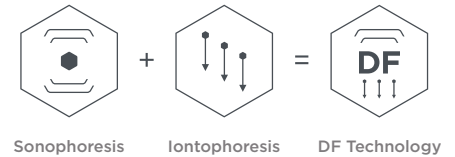

Performance Intelligence is the ultimate in skincare innovation. It combines the world's most innovative DF Technology with the power of Environ products resulting in a skin treatment system that delivers more essential nutrients to the skin where they are needed most.

