


ROLL-CIT® APPLICATION

How to use the Cosmetic Roll-CIT®/
Gold Cosmetic Roll-CIT®/
Body Cosmetic Roll-CIT®


- **STEP 1:** Start with clean skin
- **STEP 2:** Roll in a star pattern
- **STEP 3:** Apply mild pressure
- **STEP 4:**
 - **Zone 1:** Forehead and upper eyelids 1 minute
 - **Zone 2 & 3:** Cheeks and lower eyelids 1 minute each side
 - **Zone 4:** Neck and décolleté 1 minute each side
- **STEP 5:** Follow with your Environ serum of choice

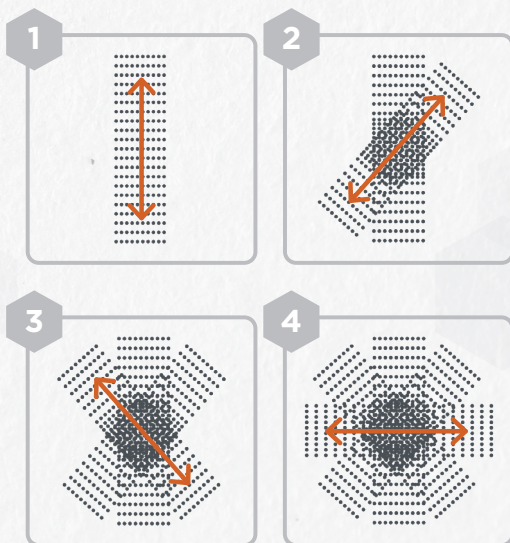

FOCUS-CIT® APPLICATION

How to use the Focus-CIT®


- **STEP 1:** Start with clean skin
- **STEP 2:** Stamp in overlapping positions
- **STEP 3:** Apply mild pressure
- **STEP 4:**
 - **Zone 1:** Forehead and upper eyelids 1 minute
 - **Zone 2:** Cheeks and lower eyelids 1 minute each side
- **STEP 5:** Follow with your Environ serum of choice

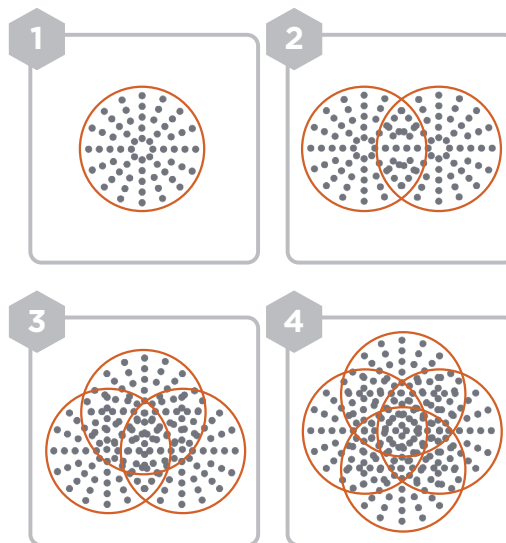

SKIN TECH+ *Range*®

WHAT DOES THE SKIN TECH+
RANGE® DO?

Environ instruments can be introduced into the evening skincare routine, and are most effective when used in combination with various Environ serums – formulated with highly intelligent ingredient combinations for the visible signs of:

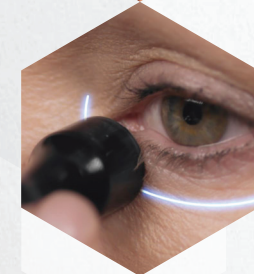

1. **Fine lines and wrinkles:** Retinol Serum, Avance Elixir, Focus Frown
2. **Dehydration:** Hydrating Serum
3. **Pigmentation:** Mela-Fade Serum System
4. **Scarring:** Retinol Serum
5. **Sensitivity:** Colostrum Gel

WHO SHOULD USE THE SKIN TECH+ RANGE®

Anyone who would use the Skin Tech+ Range is willing to invest in a skincare program that will enhance the results and skin benefits of Environ Skin Care products.

The key benefits of the Skin Tech+ Range include:

- ⬡ **Safe and effective:** Delivers more essential nutrients and intelligent ingredient combinations to where the skin needs them most
- ⬡ **Simple and easy to use:** Meticulously crafted medical-grade quality cosmetic micro-needling instruments
- ⬡ **Convenient:** Use at home and in professional treatments
- ⬡ **Minimally invasive treatments:** Maintains the benefits of the skin's integrity
- ⬡ **Scientifically proven:** Maintains the effects of the skin's integrity to achieve the healthiest-looking skin possible

BENEFITS OF THE SKIN TECH+ RANGE®

Environ's philosophy on Electro-Sonic and Micro-Needling Technology remains the same — it can be introduced only when the skin is comfortable with medium to high levels of vitamin A.

PIONEERING TECHNOLOGY

MEDICAL-GRADE QUALITY

THE ROUTE TO THE BEST RESULTS IN THE WORLD



Instruments

The Skin Tech+ Range consists of the following products:

- ⬡ Electro-Sonic DF Mobile
- ⬡ Gold Cosmetic Roll-CIT®
- ⬡ Body Cosmetic Roll-CIT®
- ⬡ Cosmetic Roll-CIT®
- ⬡ Cosmetic Focus-CIT®
- ⬡ Instrument Cleaning Solution



www.environskincare.com

ENVIRON®
a beautiful skin for a lifetime