INSTRUCTIONS FOR USE

Electro-Sonic

Ionzyme[®]
DF Machine

The Electro-Sonic Ionzyme[®] DF Machine uses industry-leading, innovative electro-sonic technology that helps stimulate the skin, allowing Environ's powerhouse formulations to be more effective, promoting essential nutrients to the skin, where they are needed most; making a real lasting difference to skin, for life.

ENVIRON[®]
Professional

Manufactured for:

Environ Skin Care (Pty) Ltd, 14 Jan Smuts Road,
Beaconvale, Parow, Cape Town, 7500, South Africa

Made in South Africa

www.environskincare.com

22/IM/DFMIFU/02

PLEASE READ AND RETAIN THIS BOOKLET FOR FUTURE REFERENCE AND WARRANTY PURPOSES.

the **BENCHMARK** *of beautiful skin*

We believe that Skin Therapy is a Science

A Beautiful Science

And that real, lasting beauty goes deeper

We believe that as professionals

There are no quick fixes

And that we too should go deeper

With a deeper understanding of the power of vitamin A

To revive, restore and rehabilitate

A deeper conviction to use the power of Performance Intelligence

To work with skin

To deliver more of what skin needs to where skin needs it most

Because we're not just skin care professionals - we're carers

We make a real difference to skin

And to our clients' lives


A lasting difference

We make sure

Skin not only looks healthier, but is healthier

Skin not only looks beautiful, but is beautiful - for life

Because beautiful-looking skin should be for everyone - for life

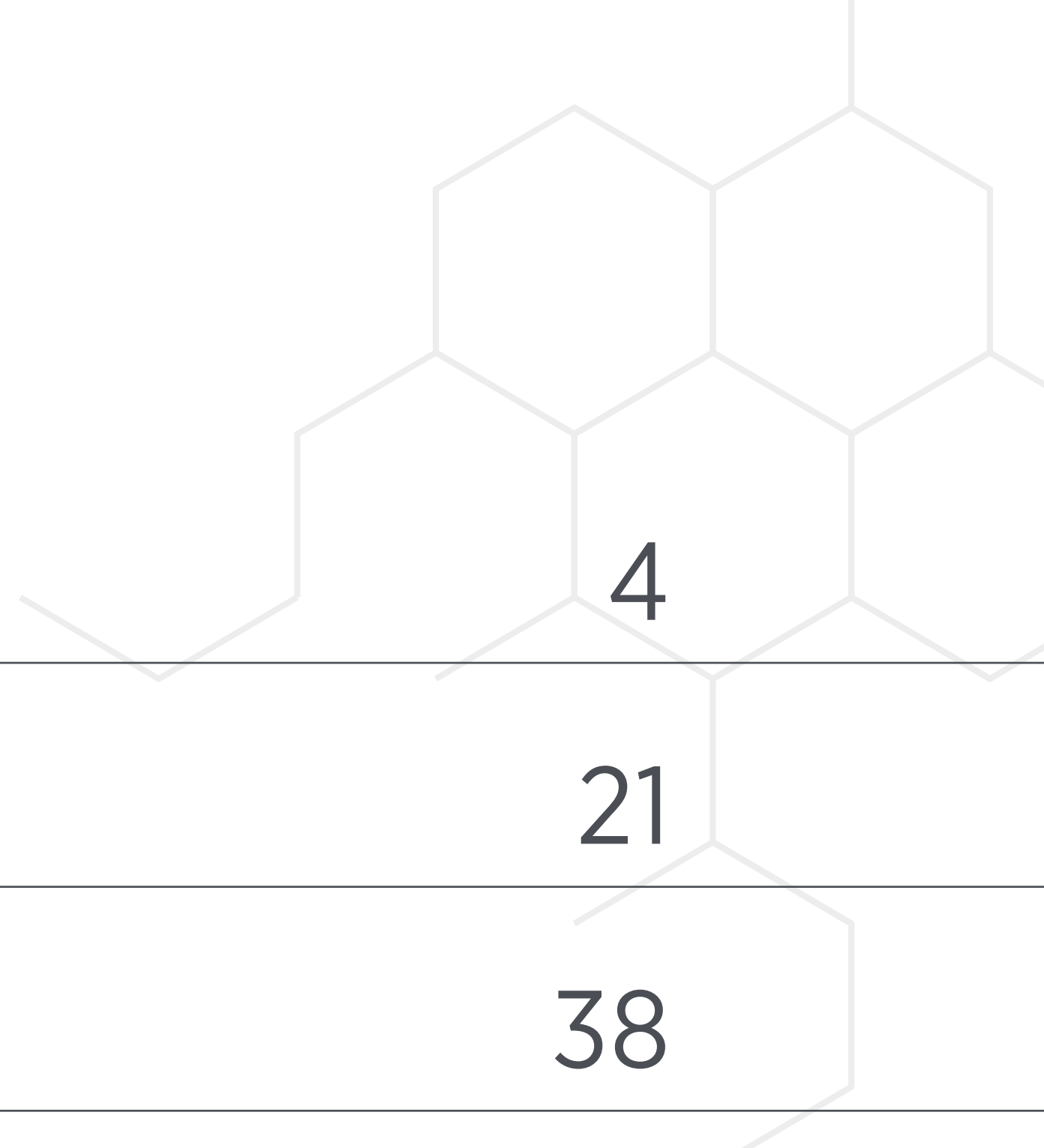ENGLISH	4
FRANÇAIS	21
DEUTSCH	38
NEDERLANDS	55
ESPAÑOL	72
PORTUGUÊS	89
ITALIANO	106
POLSKI	123
ΕΛΛΗΝΙΚΑ	140
SUOMI	157
DANSK	174
NORSK	191
한국어	208

Introduction

Congratulations on purchasing your Environ Electro-Sonic Ionzyme® DF Machine which will allow you to offer a facial like no other.



For more than 30 years, Environ has remained at the forefront of skin care innovation with revolutionary product formulations and pioneering skin care technologies, that have been proven to make a real, lasting difference to the appearance of the skin.

The revolutionary Environ Ionzyme® DF Machine, pioneered by Dr Des Fernandes, was the first in the world to combine low frequency Sonophoresis and pulsed Iontophoresis. Harnessing the power of both technologies speeds up and intensifies treatment effects. Electro-Sonic DF Technology is unique in every way and is well-regarded in the industry and can be supported by the power of an Environ vitamin A home care routine.

The Ionzyme® DF Machine is not intended for use in the diagnosis of disease or other conditions, or in the cure, mitigation, treatment or prevention of disease, and is not intended to affect the structure or any function of the living tissue and therefore is not classified as a device under the FD&C Act. 179.

Familiarize yourself with your DF Machine

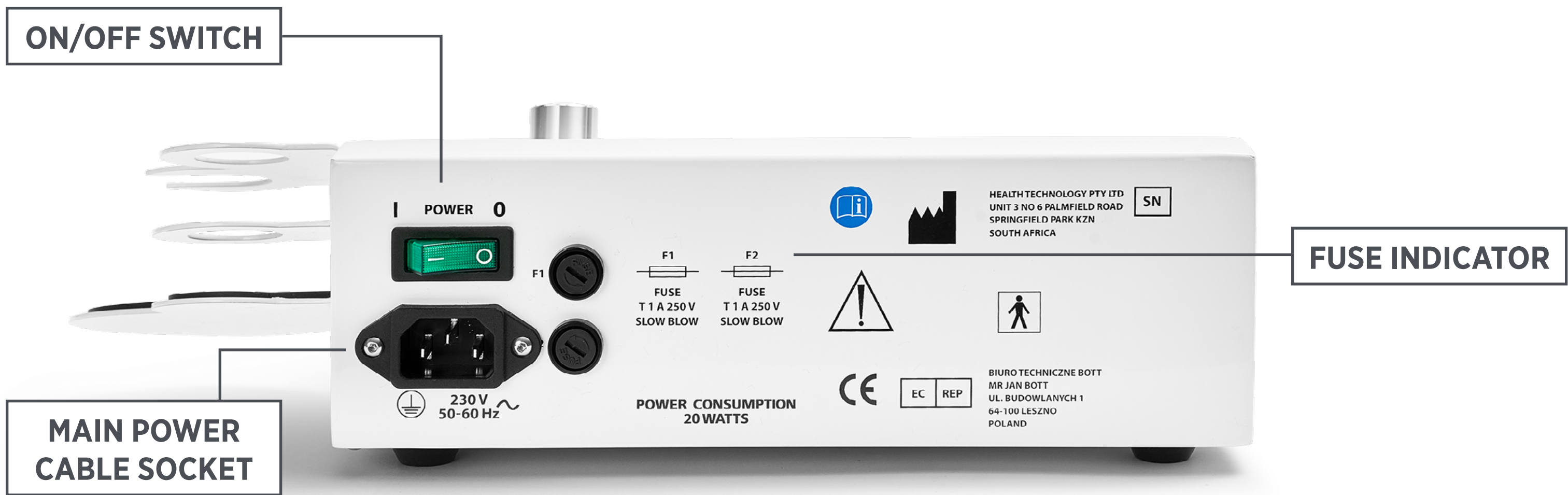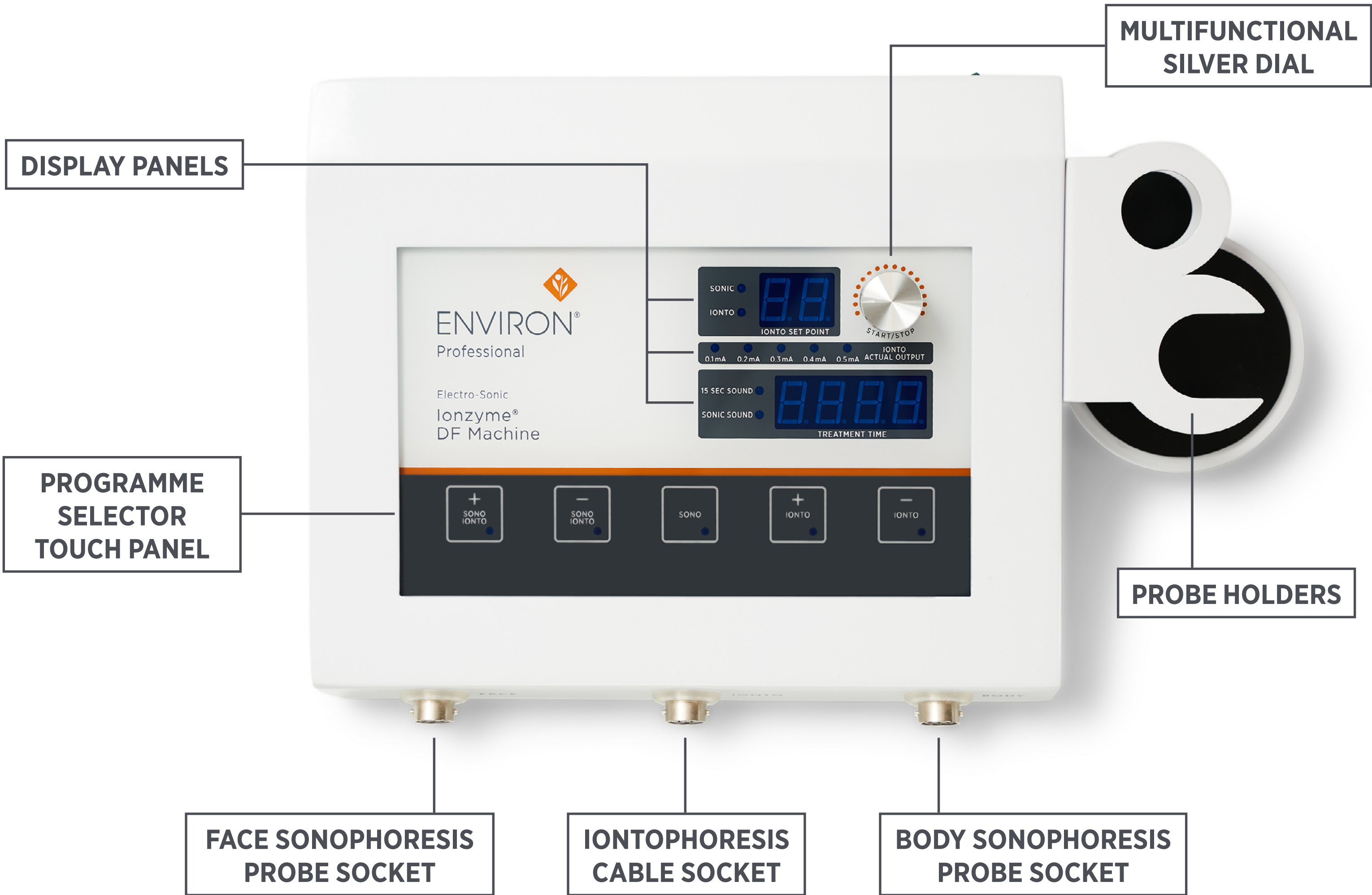
SONOPHORESIS PROBE: BODY

SONOPHORESIS PROBE: FACE

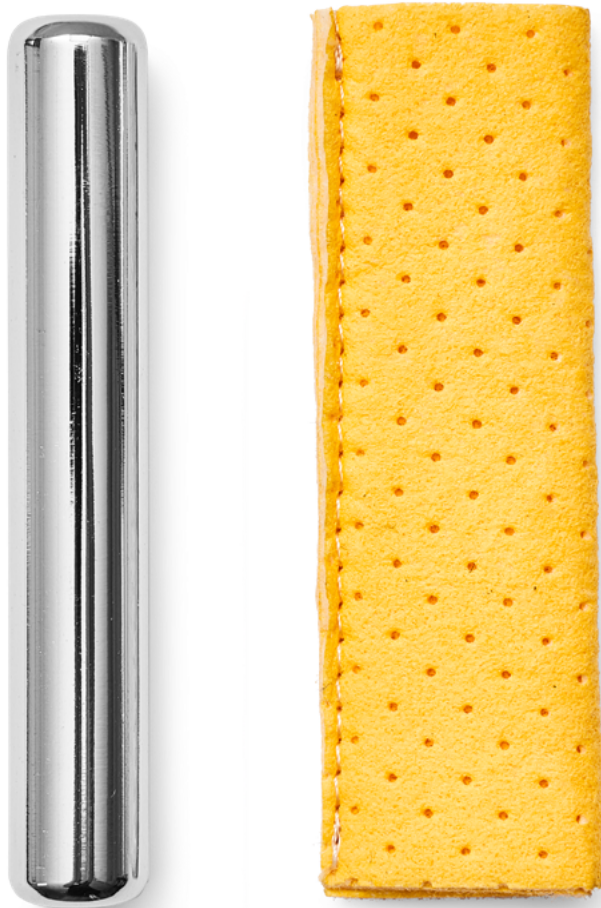**PASSIVE POLE:
ROD WITH SPONGE**


**PASSIVE POLE:
ROUND PAD AND CABLE**

**ACTIVE POLE:
CROCODILE CLIPS**


MAIN POWER CABLE

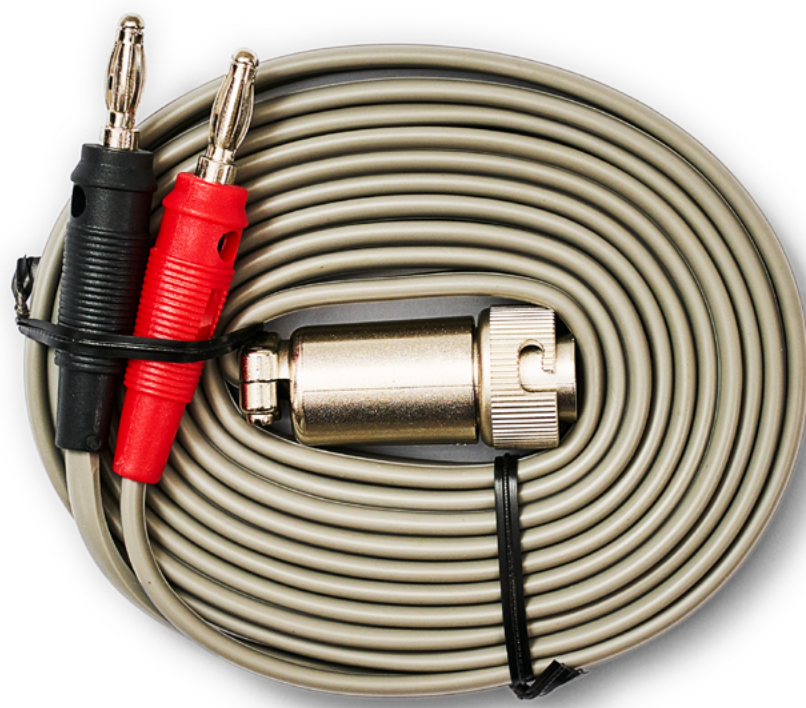IONTOPHORESIS CABLE

How does the DF Machine work

The Environ Electro-Sonic Ionzyme® DF Machine incorporates two technologies: pulsed Iontophoresis and low frequency Sonophoresis.

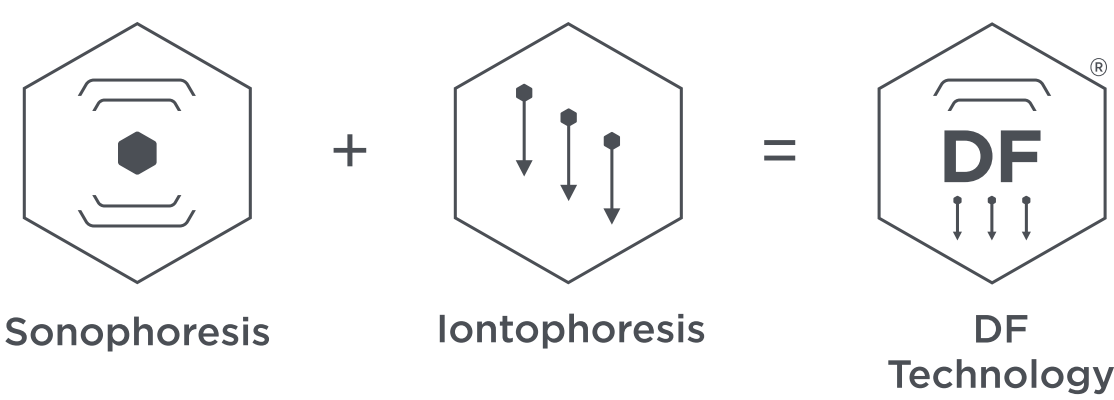Sonophoresis is a sound-based treatment system using a 22 kHz sound wave – defined as a low frequency sound wave. These sound waves are generated from the vibration of a special ceramic head and are transmitted through to the skin to disrupt the stratum corneum by means of cavitation.

Iontophoresis is a non-invasive method that makes use of a small electrical current to propel charged molecules into the skin making it more permeable for cosmetic applications. The technology is based on the principle that an electric potential will cause ions in solution to migrate according to their electrical charges.

VIEW VIDEO ON FEATURES & BENEFITS OF THE DF MACHINE

INTENDED USE

The Environ® Ionzyme® DF Machine uses low frequency Sonophoresis and Iontophoresis and is intended to temporarily disrupt the stratum corneum to improve the penetration of skincare for cosmetic applications.

The Environ® Ionzyme® DF Machine is not intended for use in the diagnosis of disease or other conditions, or in the cure, mitigation, treatment or prevention of disease, and is not intended to affect the structure or any function of the living tissue and therefore is not classified as a device under the FD&C Act. 179.

Contraindications

THIS TECHNOLOGY IS CONTRAINDICATED FOR THE FOLLOWING CONDITIONS:

1. Pacemakers and electronic implants of all kinds
2. Pregnant or lactating women
3. Epilepsy
4. Open cuts & abrasions in the treatment area
5. Cold sores (Herpes simplex type I) and ulcers.
6. Metallic implants of the face
7. Botulinum Toxin treatment and Dermal Fillers: A two week waiting period is recommended to avoid affecting their placement.
8. Sunburnt skin
9. Sufferers of angina/arrythmia

10. Haemophilia
11. Patients who have a history of stroke
12. Cancer patients or those undergoing radiotherapy

Potential Side Effects

THE FOLLOWING MIGHT BE EXPERIENCED DURING THE TREATMENT:

1. A sensation of flashing lights over the eyes during the treatment. The machine does not emit light; however, you may experience the feeling of seeing flashes of light due to the iontophoresis current.
2. A mild tingling sensation in the area that the passive pad/electrode is placed.
3. A mild tingling sensation on the skin during the treatment.
4. A high-pitched sound when the sonophoresis probe is applied to the forehead, close to ears and nasal cavities etc.
5. Excessive metallic fillings or braces may cause the client to experience a metallic taste in the mouth.

THE FOLLOWING MIGHT BE EXPERIENCED AFTER THE TREATMENT:

1. Temporary redness of the skin may be observed which subsides soon after the treatment has been completed.
2. A retinoid reaction (flaking, dryness, skin sensitivity or slight breakouts), that is due to the enhanced action of vitamin A.

No other side effects have been reported.

Warnings

1. The Electro-Sonic Ionzyme® DF Machine should only be used by a trained Environ Skin Care Professional. Environ recommends the operator of the machine should have a prior qualification in facial electricals to confidently operate the machine.
2. Use the machine only for its intended use as described in this manual. Only use the attachments recommended and provided by the manufacturer.
3. Operate the machine on the correct voltage for your country and connect this machine to a properly grounded outlet only.
4. Always attach the plug to the inlet of the machine first, then plug the cord into the wall outlet. To disconnect, turn off all controls, then remove the plug from the wall outlet.
5. The machine should never be left unattended when switched on and should always be switched off and unplugged after use. Unplug the machine from the wall outlet if it is not being used for a prolonged period. Disconnect the cord by grasping the plug. Never pull the cord itself.
6. Do not operate the machine in an area where it can be exposed to water or heat sources. Always unplug the machine if there has been any danger of water or other liquids entering the machine.

7. Never operate this machine if it has a damaged cord or plug, or if it is not working properly or if it has been dropped or damaged. Return the machine to a service centre or supplier for examination and repair.
8. Never drop or insert any object into any opening of the machine.
9. In lightning-prone areas unplug the machine from the wall outlet during an electric storm. The machine is not disconnected from the mains as long as it is plugged into the wall outlet, even if the unit itself has been turned off.
10. Close supervision is necessary when this machine is used by, on, or near children or invalids.
11. Stop the treatment immediately when irregularity or abnormality of the machine is noticed during the treatment.
12. Do not open the cabinet housing of the machine.
13. Do not place machine in a place subject to direct sunlight, excessive dust or mechanical vibration or shock.
14. Place the appliance in a location with adequate air circulation to prevent internal heat build-up in the unit.
15. Stop the treatment immediately if any pain or excessive heat is experienced.
16. Ensure that the passive pole is placed securely under the armpit, in the hand or under the lower back to ensure optimal conduction of the iontophoresis current. At all times ensure that the pole is covered with a damp cloth. Incorrect use of the passive pole will cause discomfort for the client and may result in electric burns.
17. Ensure the crocodile clips are securely attached to the conductive masque / Haize Gauze before starting the treatment. Incorrect use of the crocodile clips will cause discomfort for the client and may result in electric burns.
18. Only use Environ Skin Care products with the Electro-Sonic Ionzyme® DF Machine.

Precautions

1. Remove all jewellery before using the machine.

Warranty

1. Please retain a copy of the warranty received with the purchase of the machine. Please contact your Environ Distributor if the warranty has not been received.
2. If the DF Machine does not function according to the instructions, please discontinue use immediately and contact your Environ Distributor for assistance. Do not attempt to open and repair the machine yourself. Any interference will automatically negate the warranty.
3. The machine has a two-year warranty; the accessories have a one-year warranty.

Disposal of the DF Machine

Please return the DF Machine to your Environ Distributor for disposal who will dispose of the machine responsibly.


Setting up the DF Machine

SUPPLY VOLTAGE AND EARTHING:

The Electro-Sonic Ionzyme® DF Machine uses alternating current (AC) to power the machine. Different countries deliver different supply voltages, namely 230-volt AC and 120-volt AC. The correct machine voltage must be used in the correct country. The machine is engineered with safety in the way of an earthing system. It is important that the machine is plugged into a mains supply that has a proper earthing system. This means that the mains supply will consist of three wires, Live, Neutral and Earth. Failure to incorporate an earth will result in unstable voltages and unsafe usage for both the operator and client.

MAIN SUPPLY:

Plug in the supplied mains cable into the back of the machine. Make sure you seat the plug in as far as it will allow.

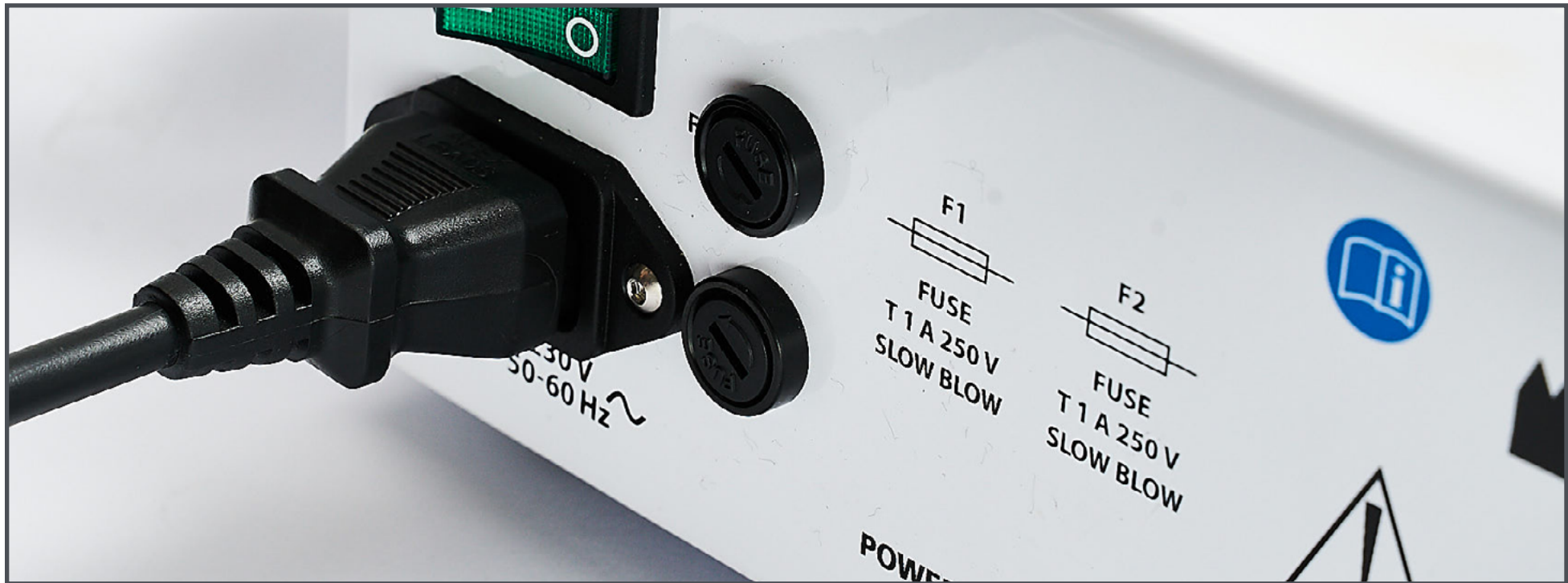FUSES:

Fuses are found at the back of the machine. These are supplied and already inserted in the machine. The fuses are there for machine protection and should not need to be changed. If the machine stops working, check the mains supply cable first. If the machine is still not getting power, then a fuse may be faulty. To check these fuses, unplug the machine from the mains and take out the fuses and check and/or replace. The fuses are specific to the machine and the correct replacement fuse ratings are indicated on the back of the machine. Fuses should only be replaced by a qualified technician. Please contact your Environ Distributor for assistance with fuse replacement.

STEP 1

Place machine on a clean, level work surface.

STEP 2

Plug the main power cable into back socket and plug into mains supply.


STEP 3

Plug in Face Sonophoresis face probe to left front socket.

STEP 4

Plug in Iontophoresis cable into centre front socket.

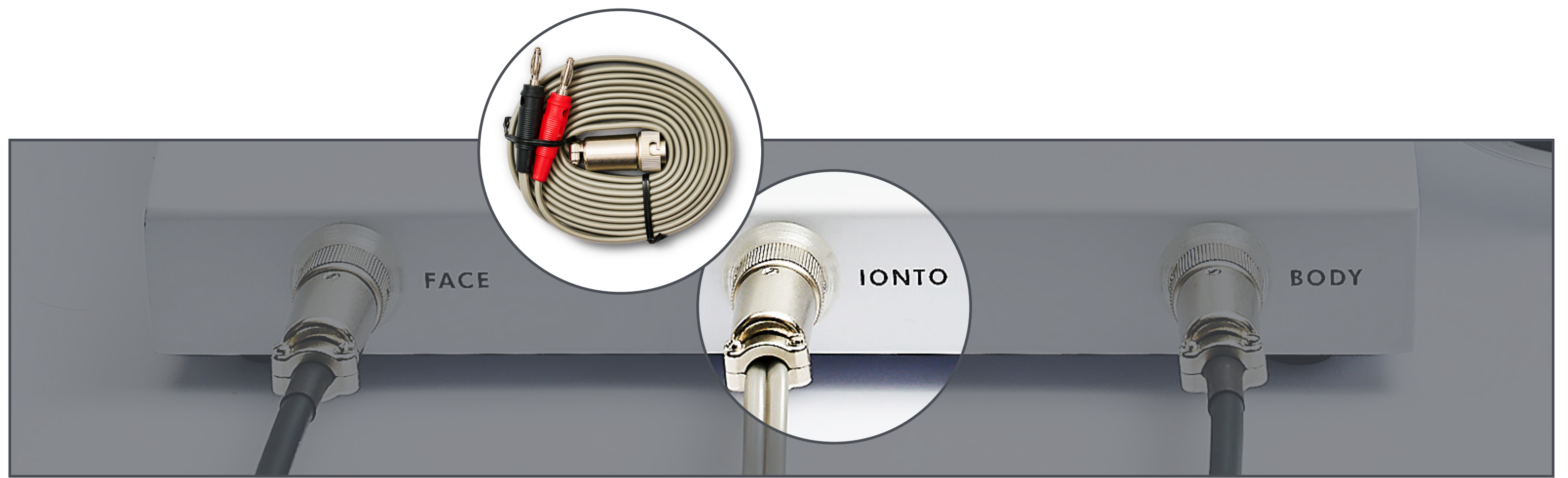STEP 5



Plug in Body Sonophoresis body probe into right front socket.

STEP 6

Plug in iontophoresis attachments into the correct iontophoresis lead outlets.

- First, ensure the iontophoresis crocodile clips are connected to the red plug.
- Next, the passive pole (pad or rod depending on preferred choice) is connected to the black plug.

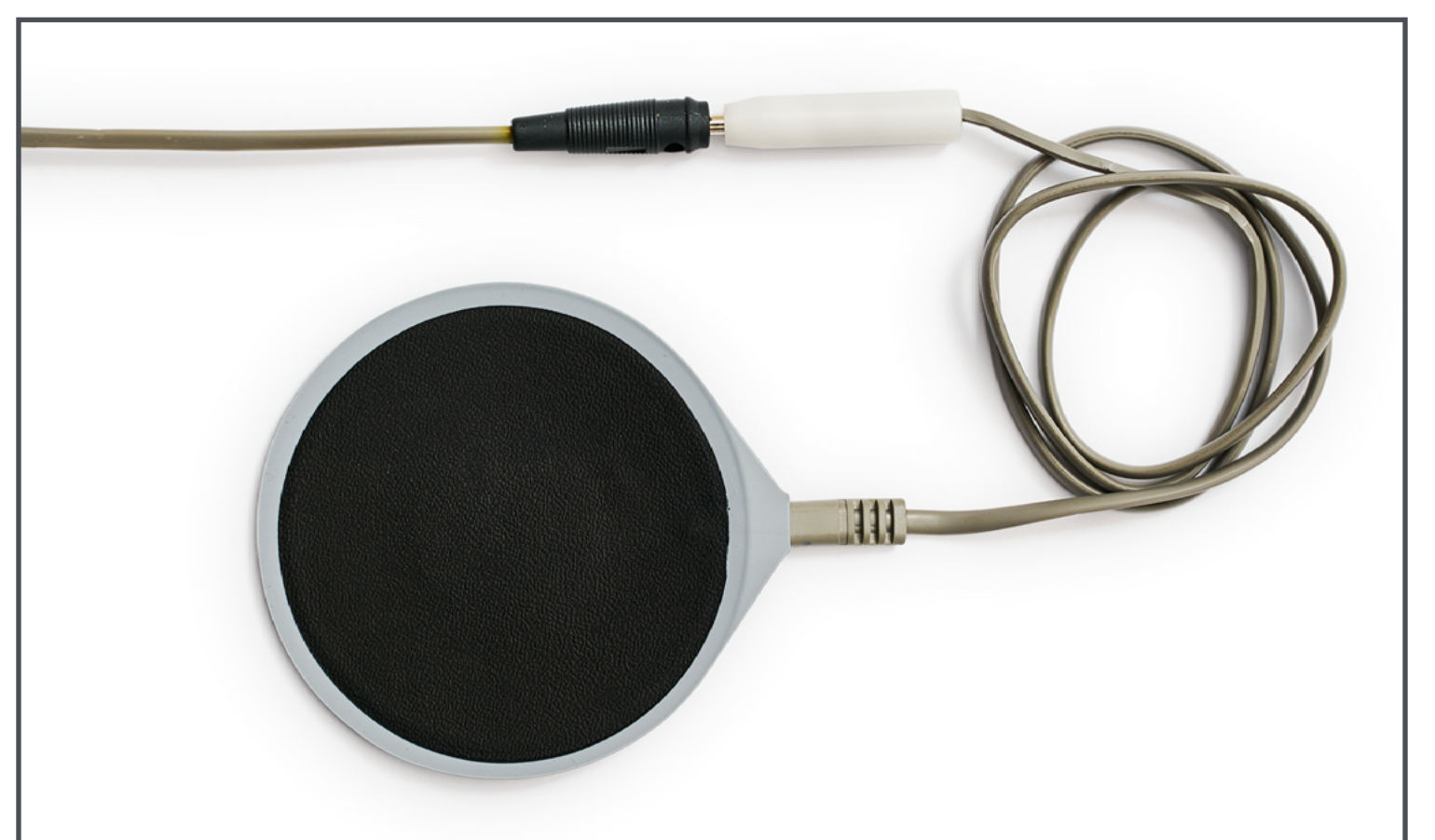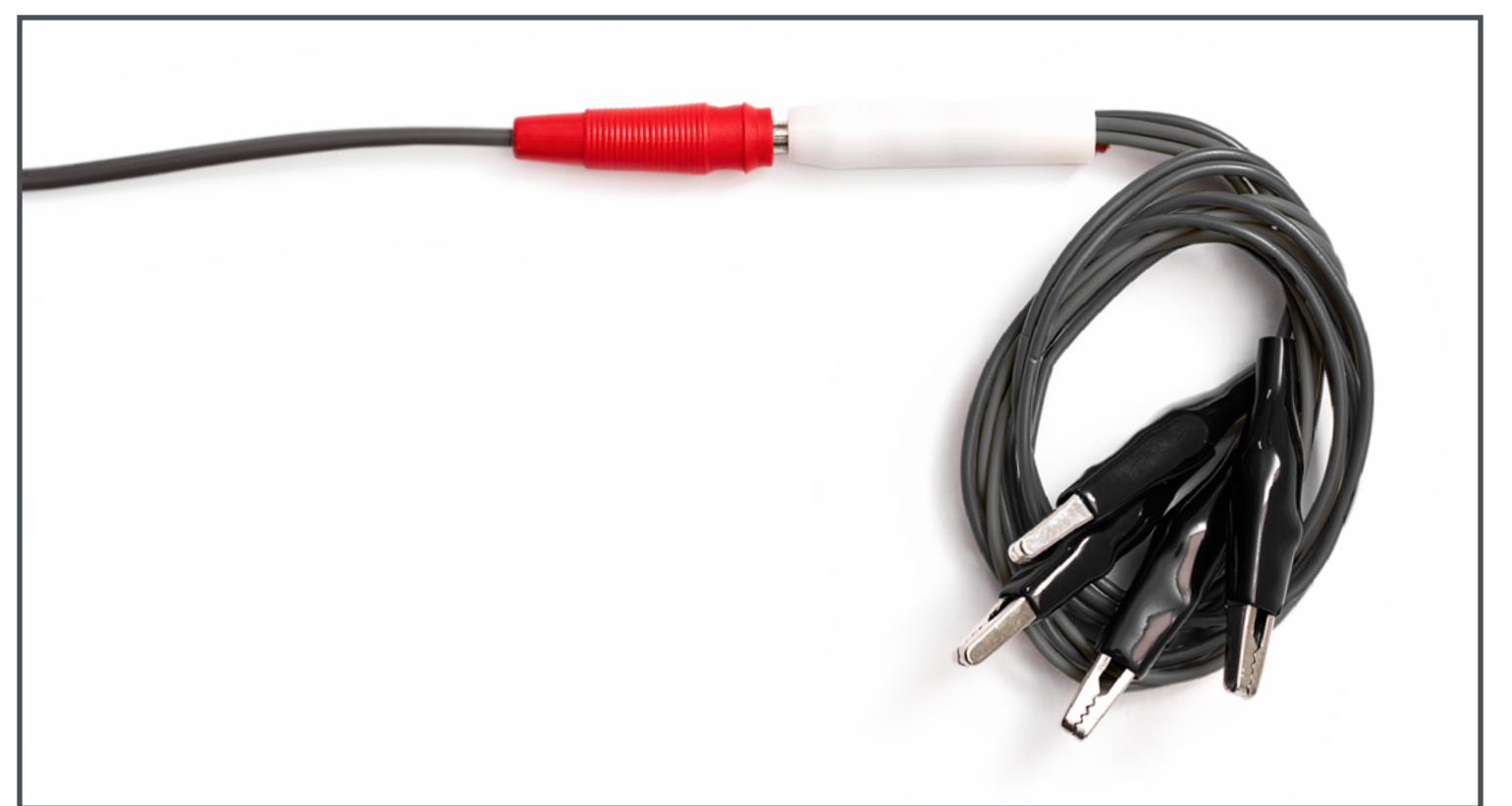Operating the DF Machine

STEP 1

SWITCH THE MACHINE ON/OFF

At the back of the machine, there is a clearly marked green mains on and off switch. This will illuminate green when the DF Machine is switched on.

At this point the DF Machine main screen will indicate default functions and settings with illuminated blue lights. These settings are set by default from the factory. Please note that the programme that is selected will also indicate at the front socket that is to be used i.e. the attachment to be used.

DEFAULT FACTORY SETTINGS

(the blue light indicates selections made)

- Programme 1: Sonophoresis and positive Iontophoresis
- 15 sec sound ON
- Sonic sound ON
- Face Sonophoresis probe output connector illuminates indicating that this will be active
- The timer will indicate 20 minutes
- Iontophoresis mA output indicator will not display until the iontophoresis current is switched on

STEP 2

PROGRAMME SETTINGS

The Electro-Sonic Ionzyme® DF Machine offers five pre-set treatment programmes:

PROGRAMME 1
Sonophoresis
and Positive
Iontophoresis
*for enhanced
penetration of
skincare*

PROGRAMME 2
Sonophoresis
and Negative
Iontophoresis
*for enhanced
cleansing*

PROGRAMME 3
Sonophoresis
*for normal
penetration of
skincare with gentle
sonic massage*


PROGRAMME 4
Positive
Iontophoresis
*for minimum
enhanced
penetration*

PROGRAMME 5
Negative
Iontophoresis
*for mild
cleansing*

CHANGING THE DEFAULT SETTINGS

To access the programming mode, press and hold the silver control dial for three seconds. All five programme lights will illuminate and flash at 1 second intervals. This indicates that it is in programming mode.

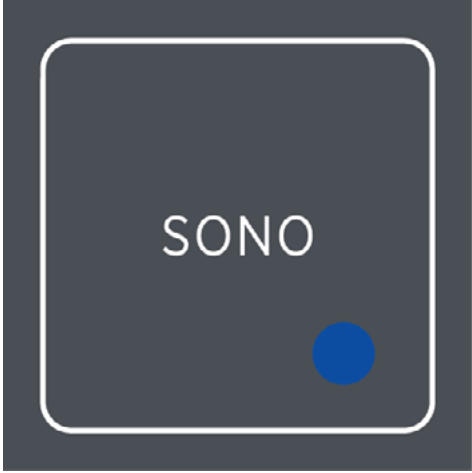INDICATION	PRESS BUTTON	RESULT
To increase the treatment time		Treatment timer increases in minutes
To decrease the treatment time		Treatment timer decreases in minutes
To turn the 15 Sec indicator sound ON or OFF		15 Sec blue indicator light will illuminate when ON and not illuminate when OFF
Turn Sonic signal sound ON or OFF		Sonic sound blue indicator light will illuminate when ON and not illuminate when OFF

After changing the settings to your desired selection press the control knob to save the new settings. These setting will now be retained as the default settings.

Please note that making setting changes could change the treatment effectiveness.

STEP 3
OPERATING THE DF MACHINE USING THE SONOPHORESIS PROBE

1. Before starting the treatment, ensure that the preferred settings are preselected.

2. Select the correct Sonophoresis treatment programme 1, 2 or 3

Please note:
If selecting programme 3, no passive pole is required as iontophoresis is not used.


3. Make sure you have placed the passive pole under the armpit, in the hand or under the lower back. Cover the passive pole with a damp cloth to ensure optimal and comfortable conduction.


4. Pressing the silver dial will start the machine and the timer will start to count down.
5. Once contact is made with the skin, the LED indicator lights will illuminate and indicate the actual current flow via the Led bar graph labelled from 0.1 to 0.5mA. Use the silver dial to increase or decrease the Iontophoresis mA output. This should be done slowly and in a controlled manner; the mA indicator will increase from 0.1 – 0.5mA in 0.1mA increments.

6. Once you have set the machine to the desired level, switch sonophoresis on by pressing the button on the probe handle. A blue indicator light will illuminate on the probe handle as well as the machine. The treatment may now start.

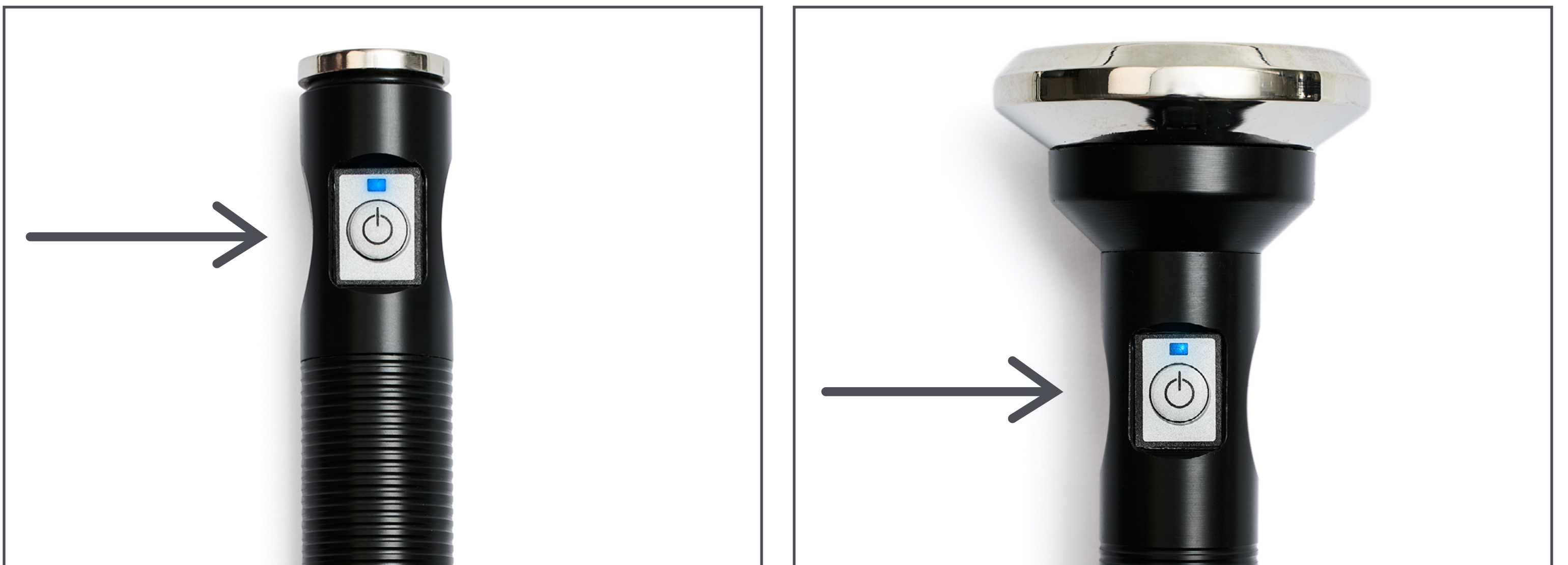- 7. Increase or decrease the Iontophoresis mA control to the desired setting at any time during the treatment by turning the silver dial either to the left (decrease) or the right (increase).
- 8. To stop the treatment, press the silver dial and the button on the sonophoresis probe. At the end of the treatment, the timer will automatically count down to 0 seconds. At this point the machine will make 5 short beeps and 'OFF' will flash on the panel. Sonophoresis and Iontophoresis will automatically stop.

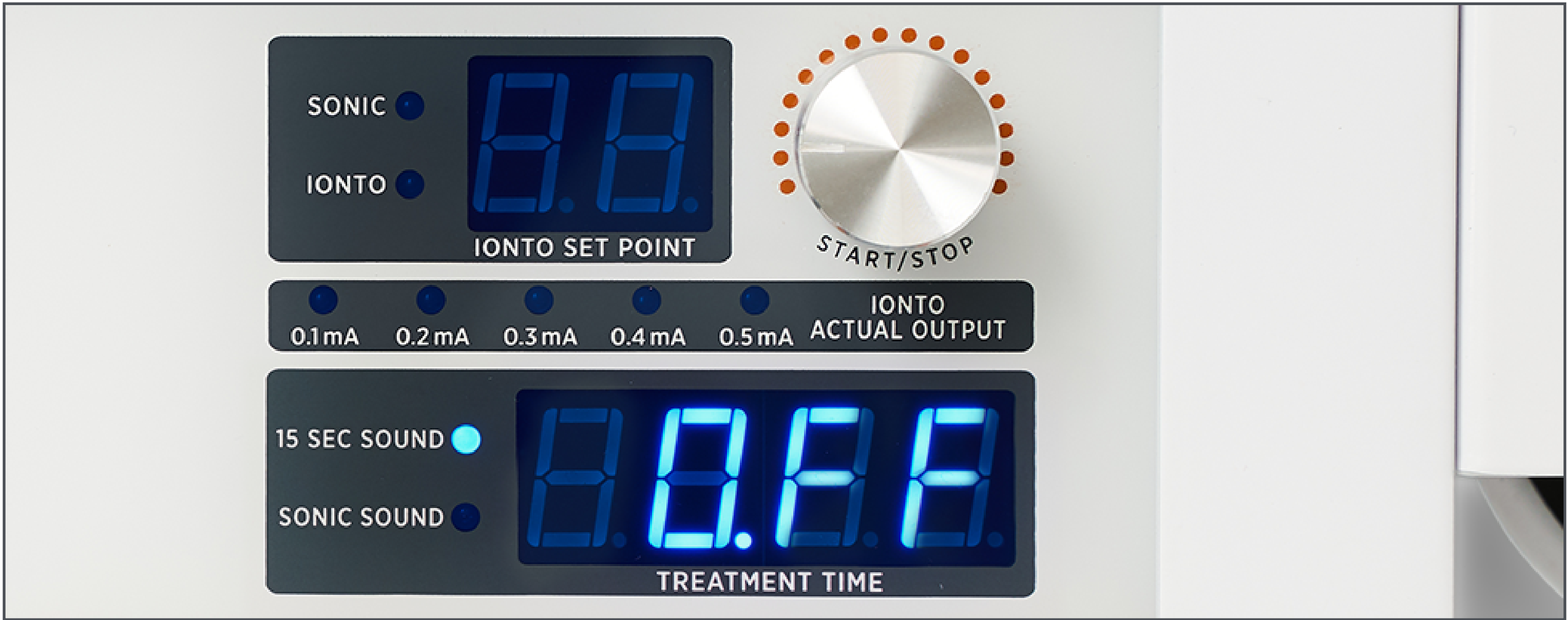
- 9. To restart the machine, simply depress the silver dial, and the machine will start again from the 20 minutes and count down in the same mode.

STEP 4
OPERATING THE DF MACHINE USING ONLY FIELD IONTOPHORESIS

- 1. Before starting the treatment, ensure that the preferred settings are preselected.


2. Select the correct Iontophoresis treatment programme 4 or 5.

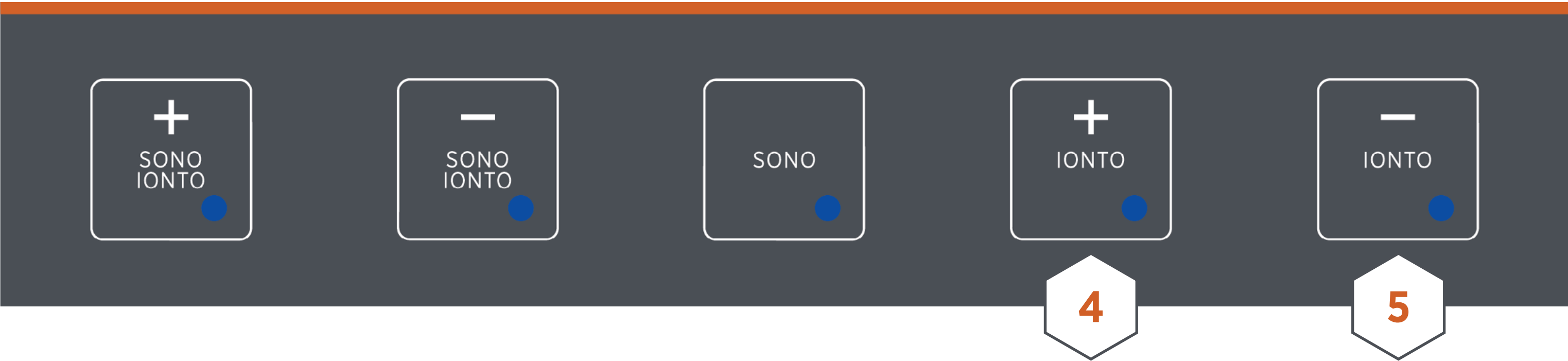
Please note:
If programme 4 or 5 are selected, the sonophoresis probe is not used in combination with the field iontophoresis, but it can be used simultaneously on a secondary area.

3. Make sure you have placed the passive pole under the armpit, in the hand or under the lower back. Cover the passive pole with a damp cloth to ensure optimal conduction.

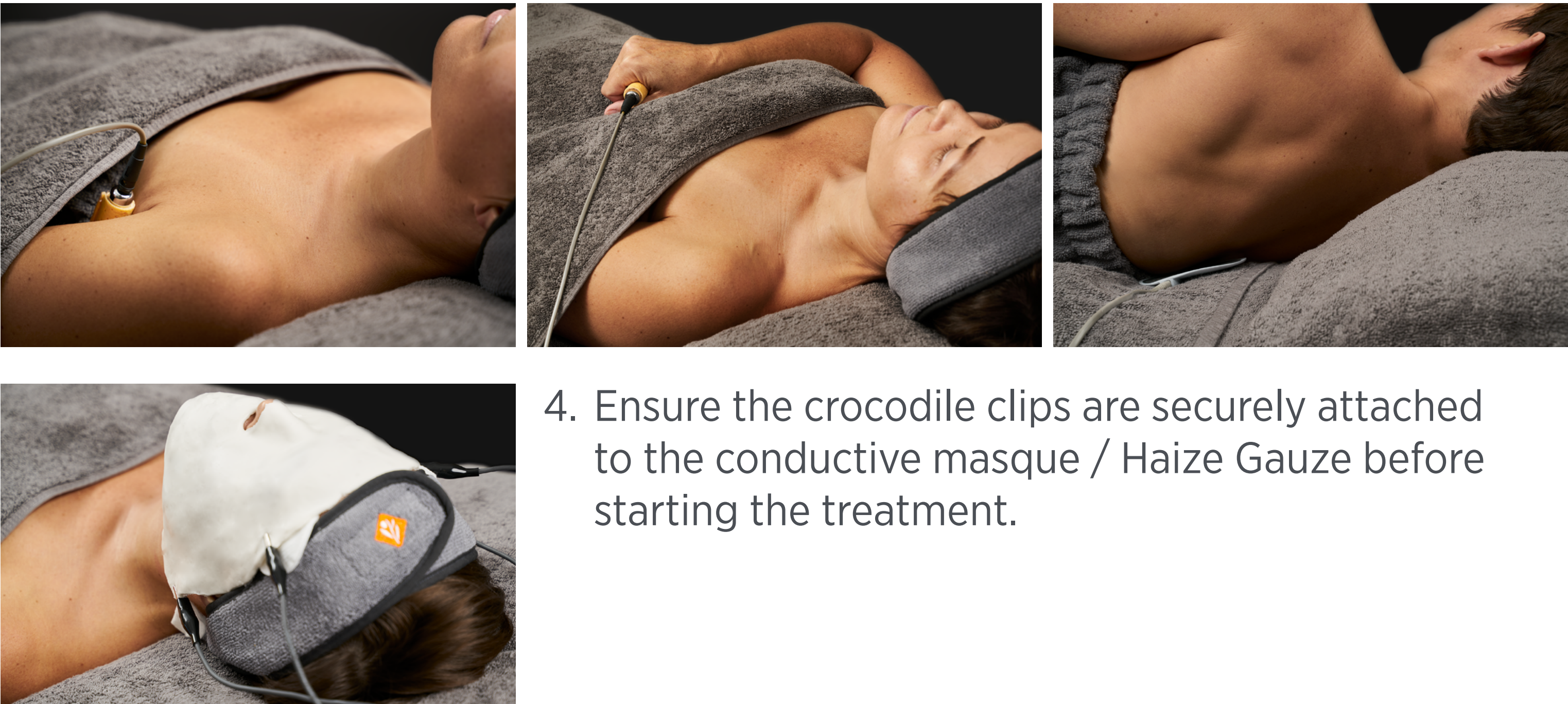4. Ensure the crocodile clips are securely attached to the conductive masque / Haize Gauze before starting the treatment.


5. Pressing the silver dial will start the machine and the timer will start to count down from 20 minutes or from what has been custom set. You may now use the silver dial to increase or decrease the Iontophoresis mA output. This should be done slowly and in a controlled manner; the mA indicator will increase from 0.1 – 0.5mA in 0.1mA increments. Once contact is made with the skin, the LED indicator lights will illuminate and indicate the actual current flow via the Led bar graph labelled from 0.1 to 0.5mA.

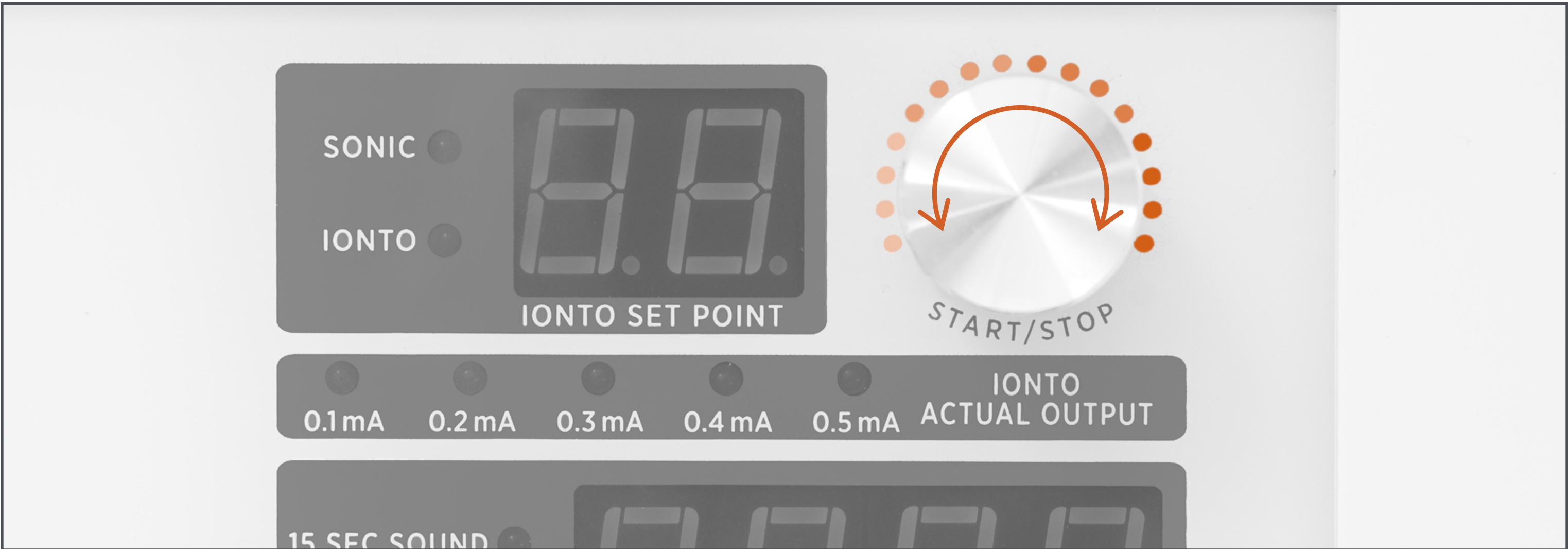6. Once you have set the machine to the desired level let the treatment continue for 20 minutes. Increase (turn to the right) or decrease (turn to the left) the Iontophoresis mA control to the desired setting at any time during the treatment. To stop the treatment, press the silver dial and the button.

7. At the end of the treatment, the timer will automatically count down to 0 seconds. At this point the DF Machine will make 5 short beeps and 'OFF' will flash on the panel. Iontophoresis will automatically stop.

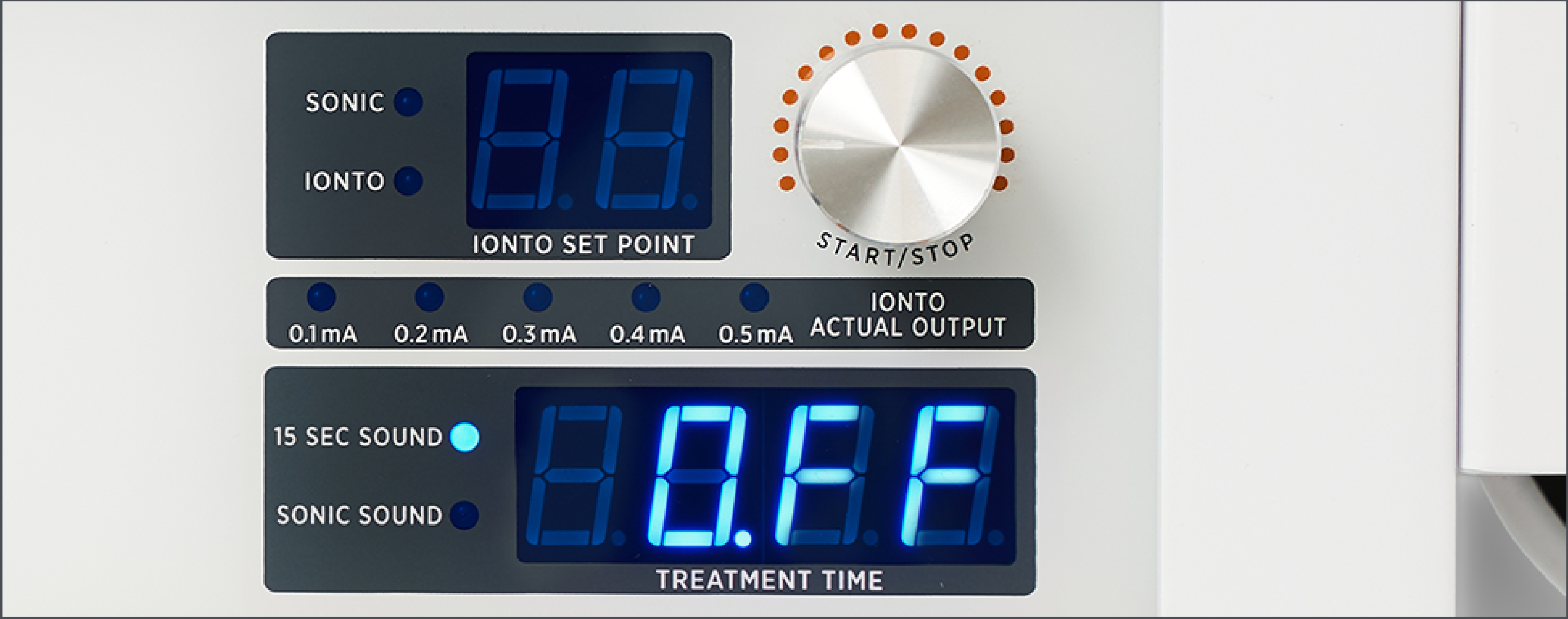8. To restart the machine, simply depress the silver dial, and the machine will start again from the 20 minutes and count down in the same mode.

*Please note:
In field Iontophoresis mode one has the option to use the Sonophoresis face or body probe at the same time as the Iontophoresis is administered but applied to a secondary area. To use either of these probes, simply pick up the required probe and depress the probe button on the handle.*

Taking care of the DF Machine

1. Gently press and / or touch the buttons on the machine.
2. To avoid scratching the surface of the machine, only use a soft, slightly damp, non-abrasive cloth to clean the machine.
3. Clean the sonophoresis probe head(s) and iontophoresis crocodile clips after every use by wiping with a soft cloth using rubbing alcohol or 75% alcohol.

Troubleshooting

PROBLEM		POSSIBLE CAUSE	SOLUTION
1.	Machine will not switch on	Incorrect or faulty connection	Ensure that the main cable is plugged into the machine and into the mains supply at the wall socket. Ensure that the power switch at the back of the machine is ON and illuminated.
2.	Machine will not switch on	Fuse	If the remedies to 1 above fail to resolve the problem, check the fuses. See Fuses section under Setting up the DF Machine. Fuses should only be replaced by a qualified technician. Please contact your Environ Distributor for assistance with fuse replacement.
3.	There is no Sonophoresis output at the face or body probe	Poor connection	Make sure that the Sonophoresis probe is plugged into the machine correctly.
4.	There is no Sonophoresis output at the face or body probe	Faulty probe	If the remedies to 3 above fail to resolve the problem and if no output is felt or heard at the probe, replace the sonophoresis probe with a new probe.
5.	There is no Iontophoresis output at the sonophoresis probe.	Poor connection	Ensure that the probe and iontophoresis lead with inactive electrode are plugged into the machine correctly. Ensure the passive pole is covered with a wet cloth and in full contact with the skin. The passive pole should be plugged into the black plug on the iontophoresis lead. Ensure you have turned on the Ion Output control to the required setting.

6.	There is no Iontophoresis at the facial clip’s outlet.	Poor connection	Ensure that the iontophoresis lead with inactive electrode is plugged into the machine correctly. Ensure the passive pole is covered with a wet cloth and in full contact with the skin. Ensure the passive pole is plugged into the black plug on the iontophoresis lead. Ensure you have turned on the Ion Output control to the required setting.
<i>Please note:</i> <i>If your machine does not function properly, please advise your supplier and do not use the machine until it has been properly serviced. All servicing and warranty work should be done by an Environ appointed service engineer.</i>			

Technical Specifications

DIMENSIONS: 340 x 230 x 90mm (Including probe bracket)
MASS: 2.9Kgs machine only
POWER: 20 Watts Max
MAXIMUM ON TIME: 30 Minutes

STORAGE AND TRANSPORT IN ORIGINAL PACKAGING

AMBIENT TEMPERATURE RANGE: - 40 °C to + 70 °C
RELATIVE HUMIDITY RANGE: 10 % to 100 %, including condensation.
ATMOSPHERIC PRESSURE RANGE: 50 kPa to 110 kPa

SONOPHORESIS

ACOUSTIC WORKING FREQUENCY: 22 kHz ±5%
EFFECTIVE RADIATING AREA FACE PROBE: 4.9 cm²
EFFECTIVE RADIATING AREA BODY PROBE: 32.1 cm²

IONTOPHORESIS

OUTPUT DC - pulsed
OUTPUT WAVE FORM - pulsed DC Square Wave with equal ON/OFF periods of 15 mS, 33Hz ±10%
OUTPUT POLARITY - either Positive or Negative (switched)
VOLTAGE - 24 V pulsed ±10% (Maximum 26.4 V)
CURRENT - 0.1mA – 0.45mA ± 10% (Maximum 0.49 mA)

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